


ZB-WAXplus Test Mix**Column:** Zebtron[™] ZB-WAXplus[™], GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** Polyethylene Glycol (PEG)**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G013-11**Oven Profile:** 140oC Isothermal**Carrier Gas:** Constant Pressure Hydrogen, 1 mL/min**Injection:** Split 100:1 1 μ L @ 250°C**Detection:** Flame Ionization (FID) (255°C)**Products used in this application:**

ZB-WAXplus Test Mix

ANALYTES:

- 1** 2-Octanone
- 2** Tetradecane
- 3** Pentadecane
- 4** 1-Octanol
- 5** Hexadecane (C16)
- 6** Methyldecanoate
- 7** Heptadecane (C17)
- 8** Methylundecanoate
- 9** 1-Decanol
- 10** Octadecane
- 11** Methyl dodecanoate
- 12** 2,6-Dimethylaniline
- 13** Nonadecane
- 14** 2,6-Dimethylphenol

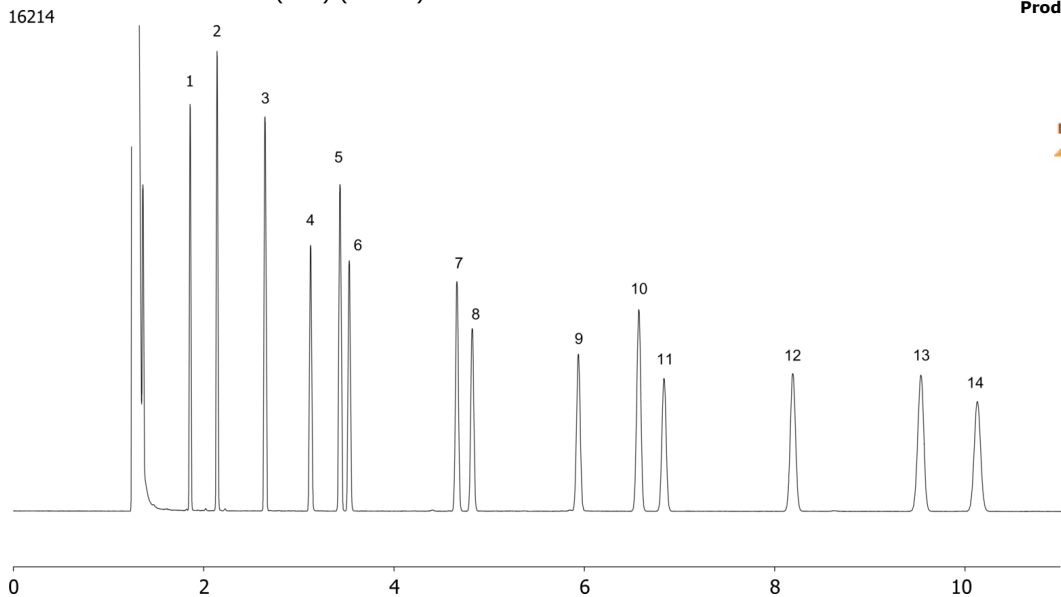