


ASTM Method D3798: Analysis of p-Xylene on ZB-WAX**Column:** Zebron[™] ZB-WAX, GC Cap. Column 30 m x 0.32 mm x 0.50 μ m, Ea**Phase:** Polyethylene Glycol (PEG)**Dimensions:** 30 meters x 0.32 mm x 0.5 μ m**Order No:** 7HM-G007-17**Oven Profile:** 75 °C to 90 °C at 3 °C/min. and hold for 5 min.**Carrier Gas:** Constant Flow Helium, 1.3 mL/min**Injection:** Split 100:1 0.1 μ L @ 220°C**Detection:** Flame Ionization (FID) (250°C)

16218

**Products used in this application:****ANALYTES:**

- 1 Toluene
- 2 Ethylbenzene
- 3 p-Xylene
- 4 m-Xylene
- 5 o-Xylene

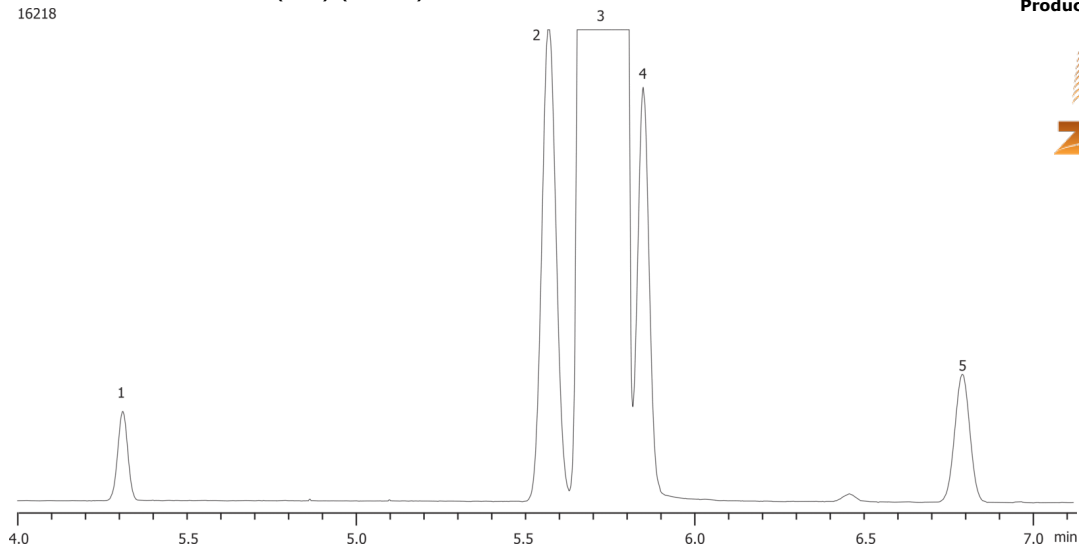