


Benzene Toluene Common Co-elutersl on ZB-1**Column:** Zebron[™] ZB-1, GC Cap. Column 100 m x 0.25 mm x 0.50 μ m, Ea**Phase:** 100% Dimethylpolysiloxane**Dimensions:** 100 meters x 0.25 mm x 0.5 μ m**Order No:** 7MG-G001-17**Oven Profile:** 60 °C to 110 °C at a ramp of 2 °C/min.**Carrier Gas:** Constant Flow Helium, 0.9 mL/min**Injection:** Split 70:1 1 μ L @ 120°C**Detection:** Flame Ionization (FID) (150°C)**Analyst Note:** 1 to 8 mg/ml Standard- 2.5 μ L 33.3 % mix per 100 μ L 2,3,3- Trimethylpentane in methanol**Products used in this application:**

16252

**ANALYTES:**

- 1** Methyl cyclopentene-1
- 2** Benzene
- 3** 2,3,3-Trimethylpentane
- 4** Toluene

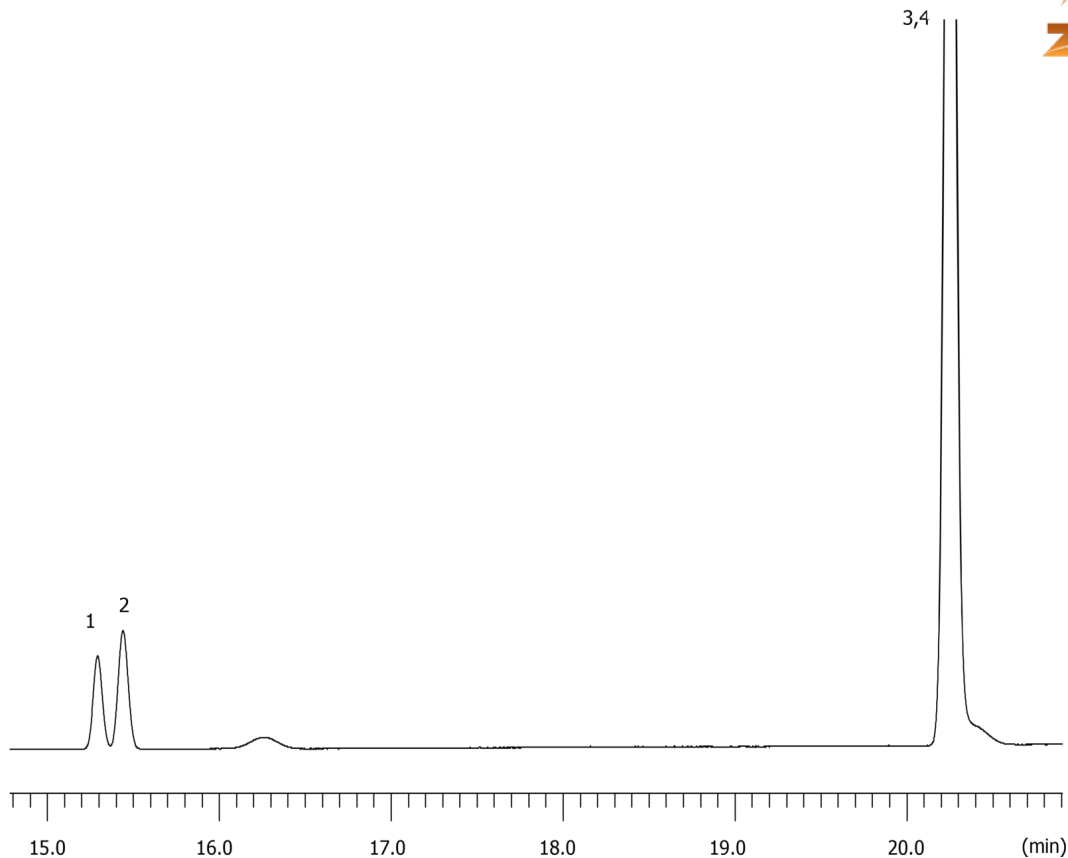