


**BTEXs on ZB-MR1****Column:** Zebtron<sup>™</sup> ZB-MultiResidue<sup>™</sup>-1, GC Cap. Column 30 m x 0.25 mm x 0.25  $\mu$ m, Ea**Phase:** Unknown or Proprietary**Dimensions:** 30 meters x 0.25 mm x 0.25  $\mu$ m**Order No:** 7HG-G016-11**Oven Profile:** 75 °C to 150 °C at 20 °C/min.**Carrier Gas:** Constant Flow Helium, 1 mL/min**Injection:** Split 200:1 0.5  $\mu$ L @ 150°C**Detection:** Flame Ionization (FID) (200°C)**Analyst Note:** Analytes are 5 mg/mL in Methanol**Products used in this application:****ANALYTES:**

- 1 Benzene
- 2 Toluene
- 3 Ethylbenzene
- 4 m-Xylene
- 5 p-Xylene
- 6 o-Xylene

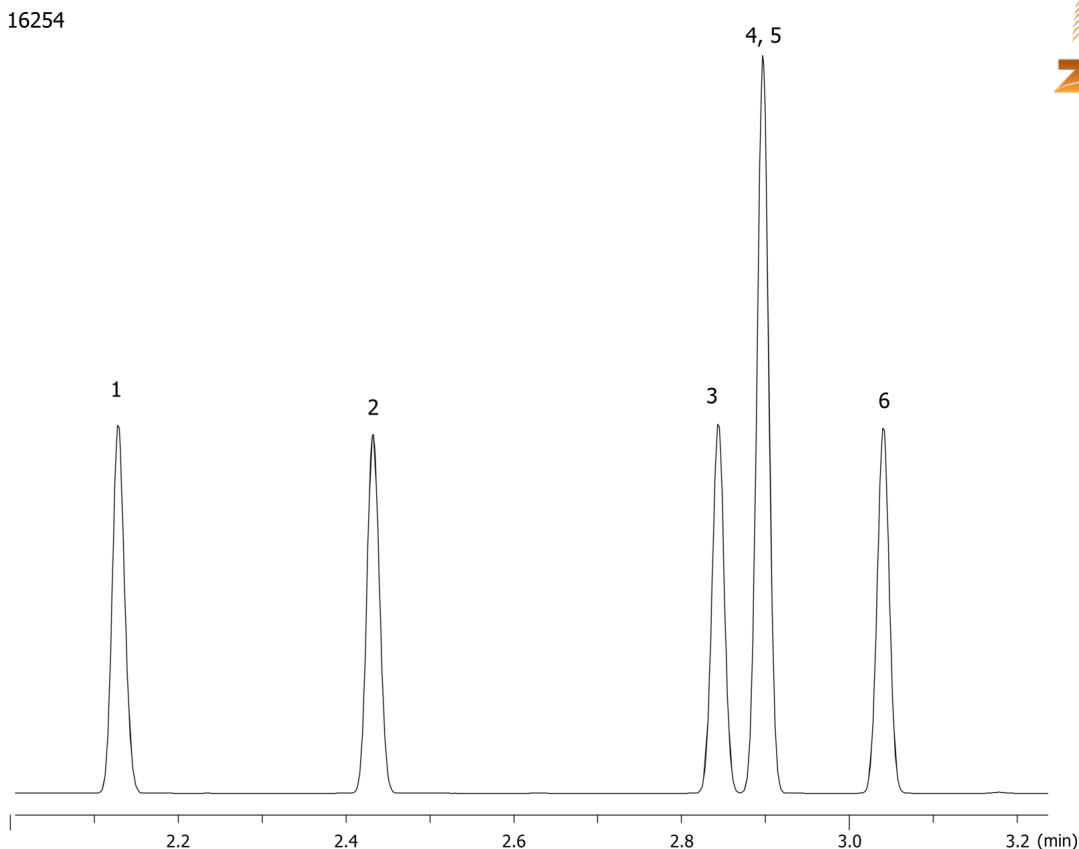