

Antihistamines on Synergi Polar-RP at low pH

Column: Synergi™ 4 µm Polar-RP 80 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4336-E0

Elution Type: Gradient

Eluent A: 20mM Potassium phosphate, pH 2.5

Eluent B: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	90	10
	2	15	20	80

Flow Rate: 1 mL/min

Col. Temp.: ambient


Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)

Detector Info: F:\Projects\0383_Nomo 5u and 3u PFP\Instrument Data\Chromatography\SW0106\PPF00314.AIA

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, POLAR-RP 4 x 3.0mm, 10/Pk Part No.: AJ0-6076

- Holder Part No.: KJ0-4282

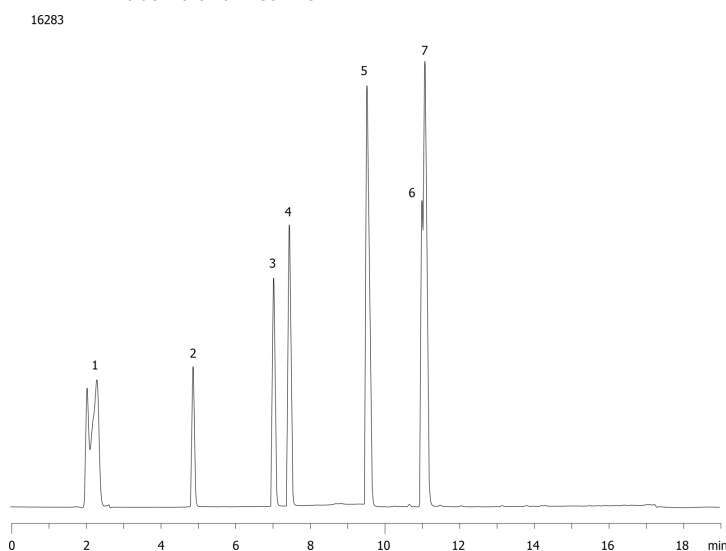

ANALYTES:

- 1 Malic acid
- 2 Doxylamine
- 3 Chlorpheniramine
- 4 Brompheniramine
- 5 Diphenhydramine
- 6 Chlorphenoxamine
- 7 Loratadine

