

Steroid Mix (11) on Luna 5um PFP(2)

Column: Luna® 5 µm PFP(2) 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4448-E0

Elution Type: Gradient

Eluent A: Water

Eluent B: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	70	30
	2	25	35	65

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, PFP 4 x 3.0 mm, 10/Pk Part No.: AJ0-8327


- Holder Part No.: KJ0-4282



Products used in this application:



16291

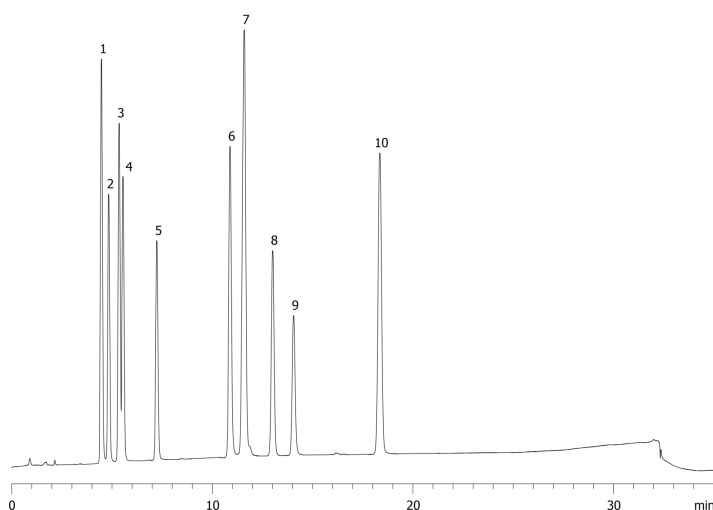

Steroid Mix (11) on Luna 5um PFP(2)

ANALYTES:

- 1** Estriol
- 2** Hydrocortisone
- 3** Prednisone
- 4** Cortisone
- 5** Corticosterone
- 6** Estradiol
- 7** Cortisone acetate
- 8** 21-Hydroxyprogesterone
- 9** 17-Hydroxyprogesterone
- 10** Estrone
- 11** Progesterone

