


Marine Oil FAMES on ZB-MR2**Column:** Zebron[™] ZB-MultiResidue[™]-2, GC Cap. Column 30 m x 0.25 mm x 0.20 μ m, Ea**Phase:** Unknown or Proprietary**Dimensions:** 30 meters x 0.25 mm x 0.2 μ m**Order No:** 7HG-G017-10**Oven Profile:** 60 °C for 2 min to 160 °C at 13 °C/min hold 5 min to 275 °C at 2 °C/min hold 2 min.**Carrier Gas:** Constant Flow Helium, 2 mL/min**Injection:** Split 5:1 1 μ L @ 250°C**Detection:** Flame Ionization (FID) (280°C)**Analyst Note:** Analytes were 200-3200ppm in methylene chloride.**Products used in this application:**

Marine Oil FAMES on ZB-MR2

ANALYTES:

- 1 Methyl myristoleate (C14:1)
- 2 Methyl myristate (C14:0)
- 3 Methyl palmitoleate (C16:1)
- 4 Methyl palmitate (C16:0)
- 5 Methyl linoleate (C18:2)
- 6 Methyl oleate (C18:1)
- 7 Methyl vaccenate (C18:1)
- 8 Methyl linolenate (C18:3)
- 9 Methyl stearate (C18:0)
- 10 Methyl arachidonate (C20:4)
- 11 Methyl eicosapentaenoate (C20:5)
- 12 Methyl 11-eicosenoate (C20:1)
- 13 Methyl 11,14-eicosadienoate (C20:2)
- 14 Methyl 11,14,17-eicosatrienoate (C20:3)
- 15 Methyl arachidate (C20:0)
- 16 Methyl docosahexaenoate (C20:6)
- 17 Methyl erucate (C22:1)
- 18 Methyl behenate (C22:0)
- 19 Methyl nervonate (C24:1)
- 20 Methyl lignocerate (C24:0)

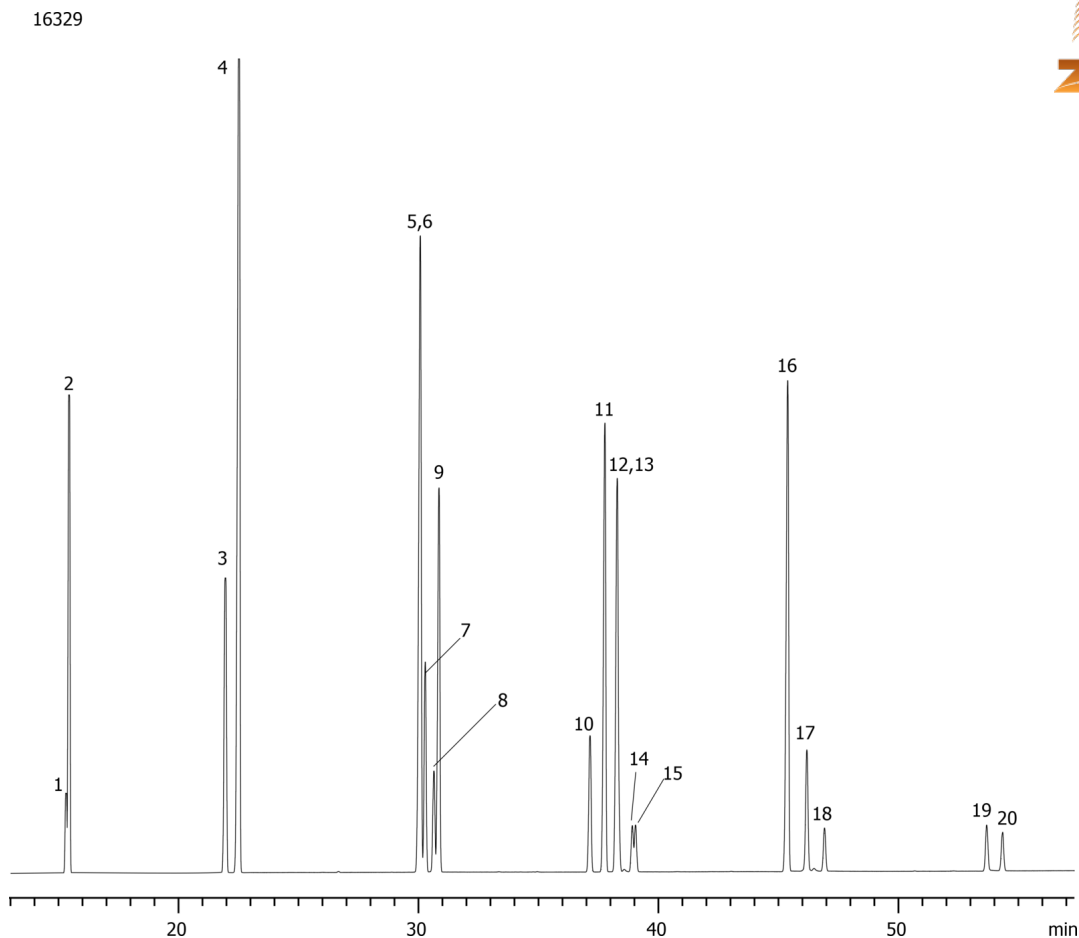