

Comparison of Water-Soluble Vitamins: Luna 5µm HILIC 80/20 ACN/Buffer, pH 4.7

Column: Luna® 5 µm HILIC 200 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4450-E0

Elution Type: Isocratic

Eluent A: Acetonitrile

Eluent B: 100mM Ammonium Acetate pH 4.7

Eluent C: 10

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B	Pct C
1	1	12	80	10	10

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Diode Array (PDA) @ 260 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.


- SecurityGuard Cartridges, HILIC 4 x 3.0 mm, 10/Pk Part No.: AJ0-8329

- Holder Part No.: KJ0-4282

16338



Products used in this application:

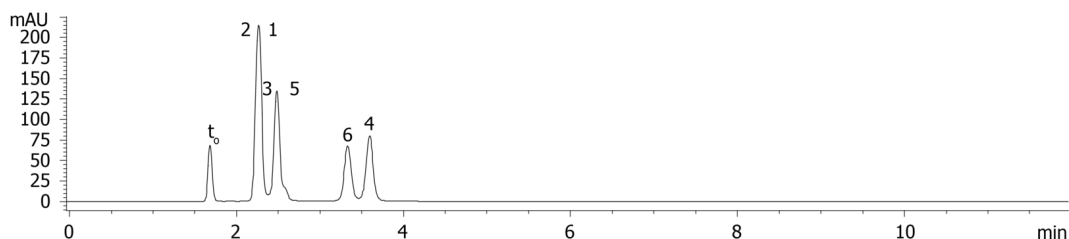

ANALYTES:

- 1 PABA
- 2 Nicotinamide
- 3 Riboflavin
- 4 Nicotinic acid
- 5 Pyridoxine
- 6 Thiamine

