


USP 467 Procedure A - Class 2 Mix A by GC/FID on ZB-624**Column:** Zebron™ ZB-624, GC Cap. Column 30 m x 0.32 mm x 1.80 µm, Ea**Phase:** 6% Cyanopropylphenyl 94% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.32 mm x 1.8 µm**Order No:** 7HM-G005-31**Oven Profile:** 40 °C for 20 min to 240 °C @ 10 °C/min for 20 min**Carrier Gas:** Constant Flow Helium, 35 cm/sec**Injection:** Headspace 5:1 1 µL @ 140°C**Detection:** Flame Ionization (FID) (250°C)**Analyst Note:** Sample Prepared as per USP Guidelines**Products used in this application:**

USP 467 Procedure A - Class 2 Mix A by GC/FID on ZB-624

ANALYTES:

- 1 Methanol
- 2 Acetonitrile
- 3 Dichloromethane
- 4 cis-1,2-Dichloroethene
- 5 Tetrahydrofuran
- 6 Cyclohexane
- 7 Methylcyclohexane
- 8 1,4-Dioxane
- 9 Toluene
- 10 Chlorobenzene
- 11 Ethylbenzene
- 12 m-Xylene
- 13 p-Xylene
- 14 o-Xylene

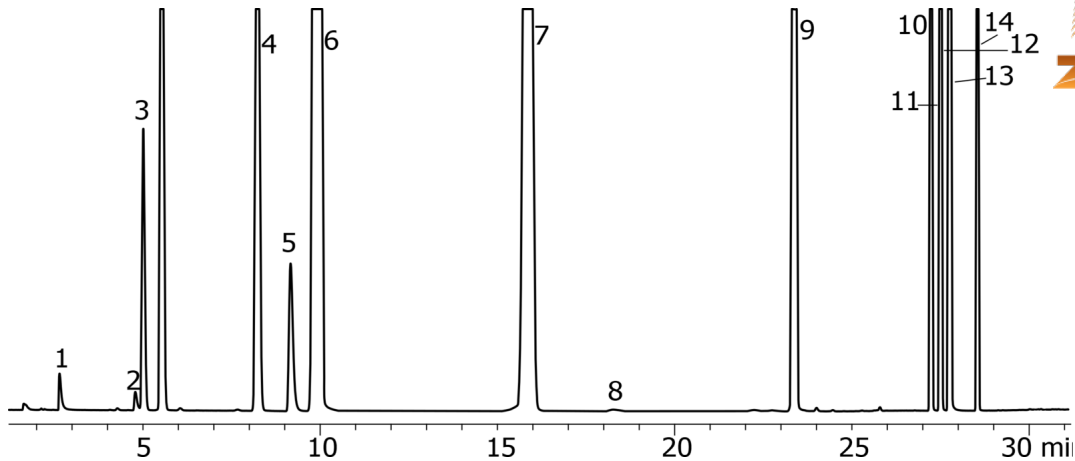