


USP 467 Procedure B - Class 2 Mix A by GC/FID on ZB-WAXplus**Column:** ZebtronTM ZB-WAXplusTM, GC Cap. Column 30 m x 0.32 mm x 0.25 μ m, Ea**Phase:** Polyethylene Glycol (PEG)**Dimensions:** 30 meters x 0.32 mm x 0.25 μ m**Order No:** 7HM-G013-11**Oven Profile:** 50 °C for 20 min to 165 °C @ 6 °C/min for 20 min**Carrier Gas:** Constant Flow Helium, 35 cm/sec**Injection:** Headspace 5:1 1 μ L @ 140°C**Detection:** Flame Ionization (FID) (250°C)**Analyst Note:** Sample prepared as per USP Method**Products used in this application:**

USP 467 Procedure B - Class 2 Mix A by GC/FID on ZB-WAXplus

ANALYTES:

- 1 Cyclohexane
- 2 Methylcyclohexane
- 3 Tetrahydrofuran
- 4 Methanol
- 5 Dichloromethane
- 6 cis-1,2-Dichloroethene
- 7 Acetonitrile
- 8 Toluene
- 9 1,4-Dioxane
- 10 Ethylbenzene
- 11 p-Xylene
- 12 m-Xylene
- 13 o-Xylene
- 14 Chlorobenzene

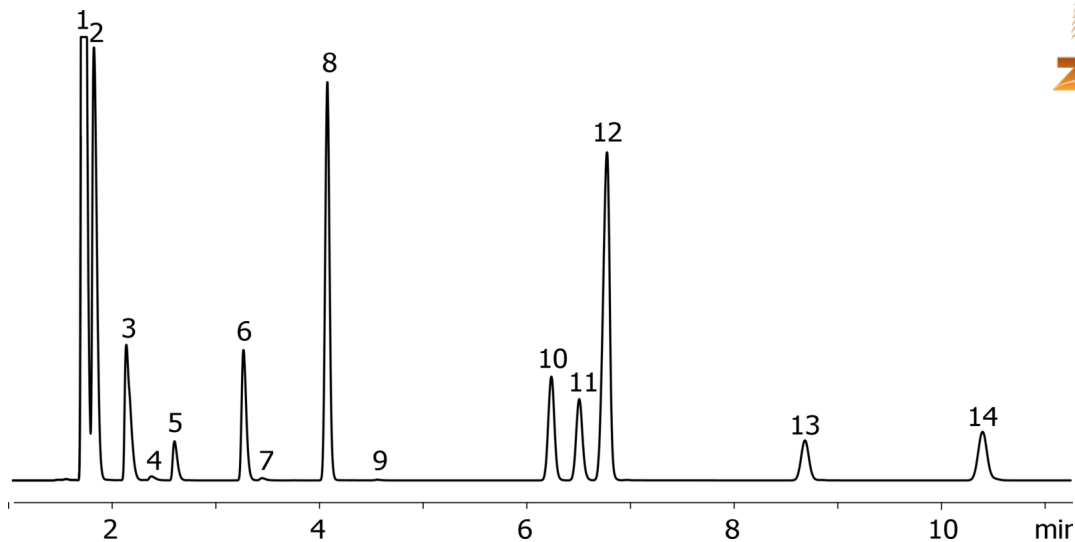