


USP 467 Procedure B - Class 2 Mix B by GC/FID on ZB-WAXplus**Column:** Zebtron[™] ZB-WAXplus[™], GC Cap. Column 30 m x 0.32 mm x 0.25 μ m, Ea**Phase:** Polyethylene Glycol (PEG)**Dimensions:** 30 meters x 0.32 mm x 0.25 μ m**Order No:** 7HM-G013-11**Oven Profile:** 50 °C for 20 min to 165 °C @ 6 °C/min for 20 min**Carrier Gas:** Constant Flow Helium, 35 cm/sec**Injection:** Headspace 5:1 1 μ L @ 140°C**Detection:** Flame Ionization (FID) (250°C)**Analyst Note:** Sample prepared as per USP Method**Products used in this application:****ANALYTES:**

- 1 Hexane
- 2 1,2-Dimethoxyethane
- 3 Trichloroethylene
- 4 Chloroform
- 5 Methyl butyl ketone
- 6 Nitromethane
- 7 Pyridine
- 8 Tetralin

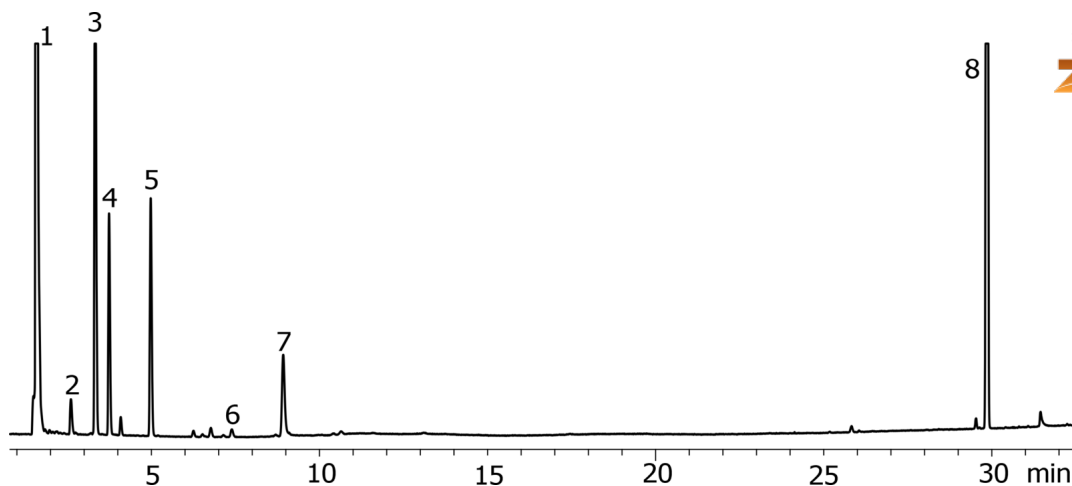