


Cognacs: Courvoisier® on ZB-WAXplus**Column:** Zebron™ ZB-WAXplus™, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** Polyethylene Glycol (PEG)**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G013-11**Oven Profile:** 35 °C for 7 min t to 60 °C @ 5 °C/min for 2 min to 210 °C @ 10 °C.**Carrier Gas:** Constant Flow Hydrogen, 1 mL/min**Injection:** Headspace 25:1 200 µL @ 210°C**Detection:** Flame Ionization (FID) (230°C)**Analyst Note:** Static Headsapce injection of 1 ml Courvoisier in 20 ml headspace sample vial at 80 °C for 20 min.**Products used in this application:**

16358

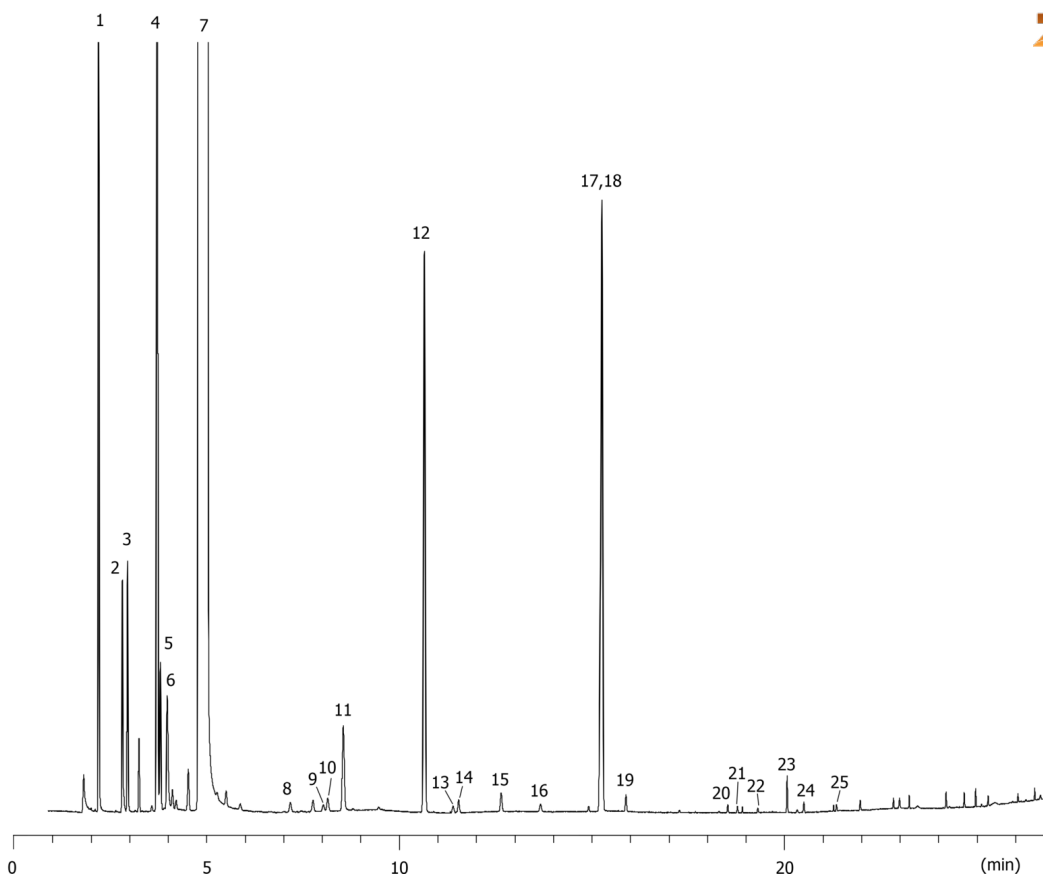

Cognacs: Courvoisier® on ZB-WAXplus

ANALYTES:

- 1 Acetaldehyde
- 2 Isobutanal
- 3 Ethyl formate
- 4 Ethyl acetate
- 5 Acetal
- 6 Methanol
- 7 Ethanol
- 8 Isobutyl acetate
- 9 2-Butanol
- 10 Ethyl butanoate
- 11 1-Propanol
- 12 Isobutanol
- 13 Allyl alcohol
- 14 Isoamyl acetate
- 15 1-Butanol
- 16 4-Methyl-2-pentanol
- 17 Methyl-2-butanol
- 18 Methyl-3-butanol
- 19 Ethyl hexanoate
- 20 Ethyl lactate
- 21 Hexanol
- 22 cis-3-Hexenol
- 23 Ethyl octanoate
- 24 Furfural
- 25 Benzaldehyde

