


Borneol on ZB-MR1**Column:** Zebron[™] ZB-MultiResidue[™]-1, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** Unknown or Proprietary**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G016-11**Oven Profile:** 60 °C for 2 min to 270 °C @ 10 °C/min**Carrier Gas:** Constant Flow Helium, 0.3 mL/min**Injection:** Splitless :1 1 μ L @ 270°C**Detection:** Mass Selective (MSD) (230°C)**Analyst Note:** Dilute Borneol in Methanol**Products used in this application:**

16360

**ANALYTES:**

- 1 Isoborneol
- 2 Borneol

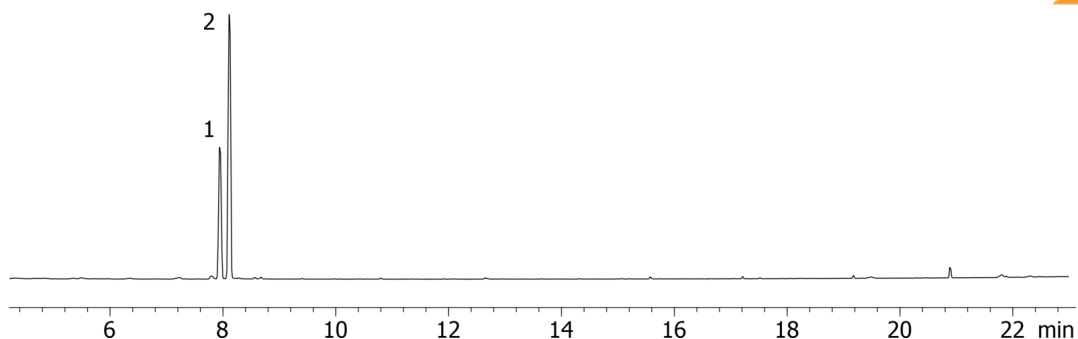