

Luna 2.5 um C18(2) 2mm ID without SecurityGuard

Column: Luna® 2.5 µm C18(2) 100 Å, LC Column 50 x 3 mm, Ea

Dimensions: 50 x 3 mm ID

Order No: 00B-4446-Y0

Elution Type: Isocratic

Eluent A: Water

Eluent B: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	35	65

Flow Rate: 0.8 mL/min

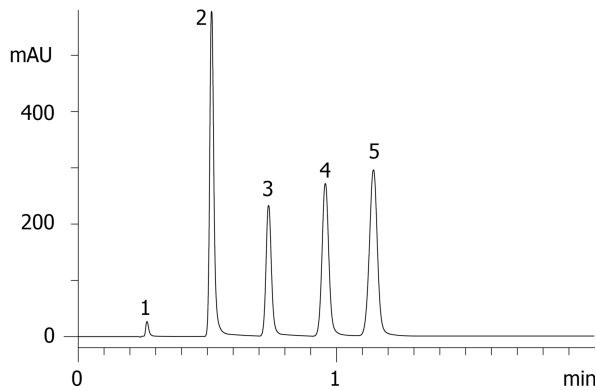
Col. Temp.: ambient

Detection: UV-Vis Abs.-Diode Array (PDA) @ (30 °C)

Detector Info: without guard column
16364



Products used in this application:

**ANALYTES:**

- | | | |
|---|--------------|---------------------------|
| 1 | Uracil | Retention Time: 0.271 min |
| 2 | Acetophenone | Retention Time: 0.521 min |
| 3 | Benzene | Retention Time: 0.74 min |
| 4 | Toluene | Retention Time: 0.96 min |
| 5 | Naphthalene | Retention Time: 1.147 min |