


Nitrogen and Phosphorous Pesticides on ZB-MR1**Column:** Zebtron™ ZB-MultiResidue™-1, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** Unknown or Proprietary**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G016-11**Oven Profile:** 60 °C for 0.5 min to 320 °C @ 10 °C/min for 5 min**Carrier Gas:** Constant Flow Helium, 0.9 mL/min**Injection:** Splitless :1 1 µL @ 260°C**Detection:** Mass Selective (MSD) (180°C)

16377

**Products used in this application:****ANALYTES:**

- 1 TEPP
- 2 Sulfotep
- 3 Monochrotophos
- 4 Dimethoate
- 5 Malathion
- 6 Parathion
- 7 EPN

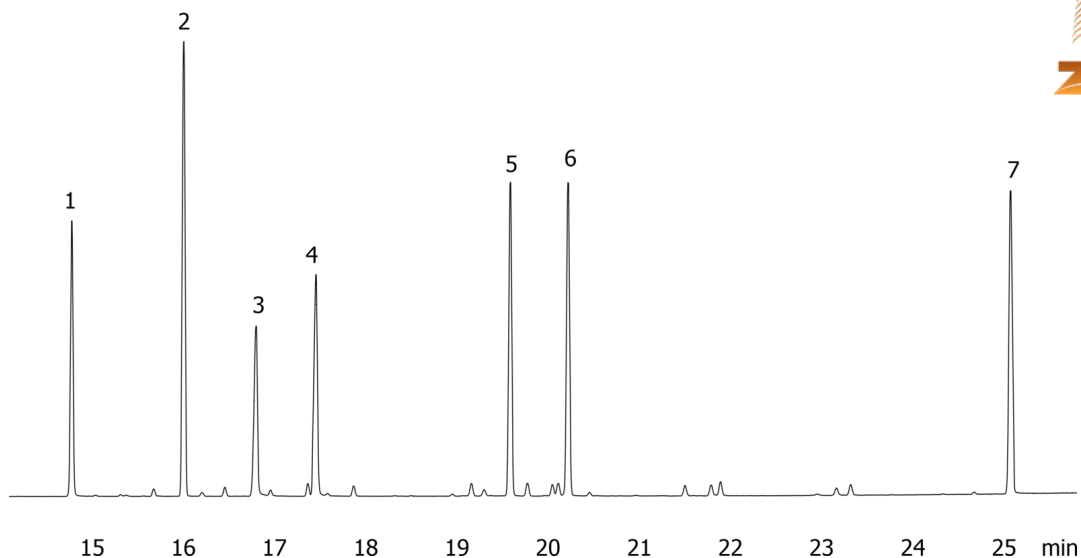