


Nitrogen and Phosphorous Pesticides on ZB-MR1**Column:** Zebtron[™] ZB-MultiResidue[™]-1, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** Unknown or Proprietary**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G016-11**Oven Profile:** 80 °C for 0.5 min to 150 °C @ 10 °C/min to 240 °C @ 4 °C/min to 320 °C @ 15 °C/min for 3 min**Carrier Gas:** Constant Flow Helium, 0.9 mL/min**Injection:** Splitless :1 1 μ L @ 260°C**Detection:** Mass Selective (MSD) (180°C)

16380

**Products used in this application:****ANALYTES:**

- 1** DEET
- 2** Terbacil
- 3** Metribuzin
- 4** Bromacil
- 5** Triadimefon
- 6** Tricyclazole
- 7** Hexazinone

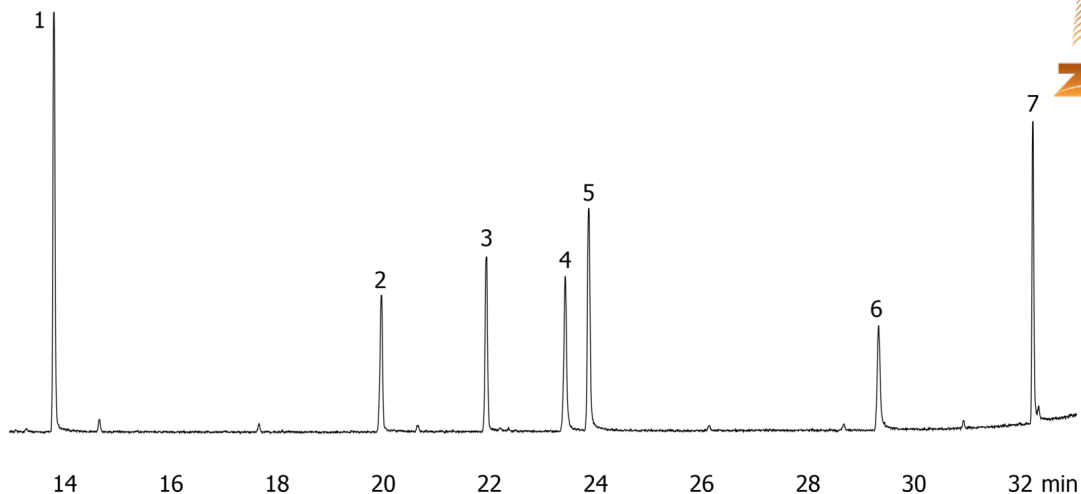