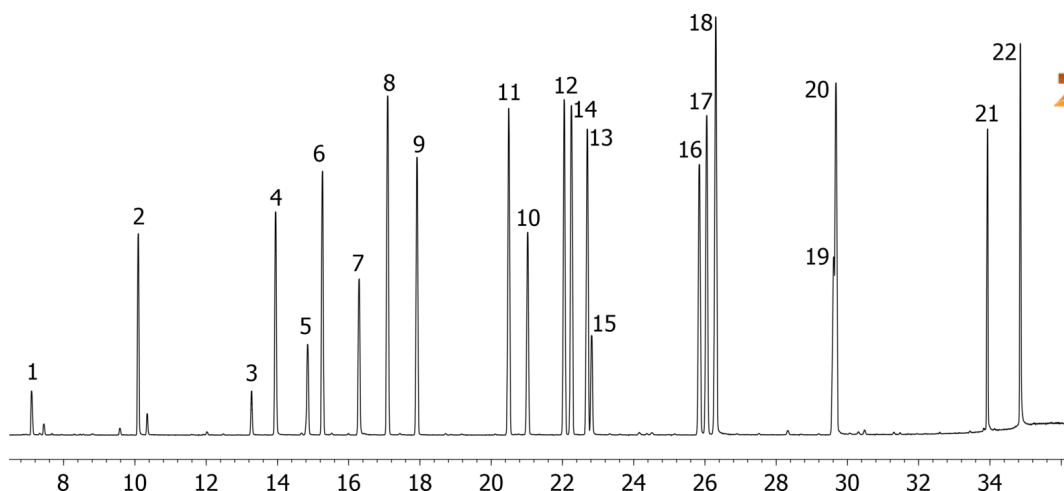


Nitrogen and Phosphorous Pesticides on ZB-MR2**Column:** Zebtron[™] ZB-MultiResidue[™]-2, GC Cap. Column 30 m x 0.25 mm x 0.20 µm, Ea**Phase:** Unknown or Proprietary**Dimensions:** 30 meters x 0.25 mm x 0.2 µm**Order No:** 7HG-G017-10**Oven Profile:** 80 °C for 0.5 min to 150 °C @ 10 °C/min to 240 °C @ 4 °C/min to 320 °C @ 15 °C/min for 3 min**Carrier Gas:** Constant Flow Helium, 0.9 mL/min**Injection:** Splitless :1 1 µL @ 260°C**Detection:** Mass Selective (MSD) (180°C)

16385

**Products used in this application:**

Nitrogen and Phosphorous Pesticides on ZB-MR2

ANALYTES:

- 1 Dichlorvos
- 2 Mevinphos
- 3 Demeton
- 4 Ethoprop
- 5 Naled
- 6 Phorate
- 7 Demeton
- 8 Diazinon
- 9 Disulfoton
- 10 Methyl parathion
- 11 Ronnel
- 12 Chlorpyrifos
- 13 Fenthion
- 14 Trichloronate
- 15 Merphos
- 16 Tetrachlorvinphos
- 17 Tokuthion
- 18 Merphos oxone
- 19 Fensulfothion
- 20 Bolstar
- 21 Azinphos-methyl
- 22 Coumaphos

