

Oligo WAX - AP fragment & DNA 20mer**Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** Water**Eluent B:** Acetonitrile**Eluent C:** 100mM Tris pH8**Eluent D:** 2M NaCl

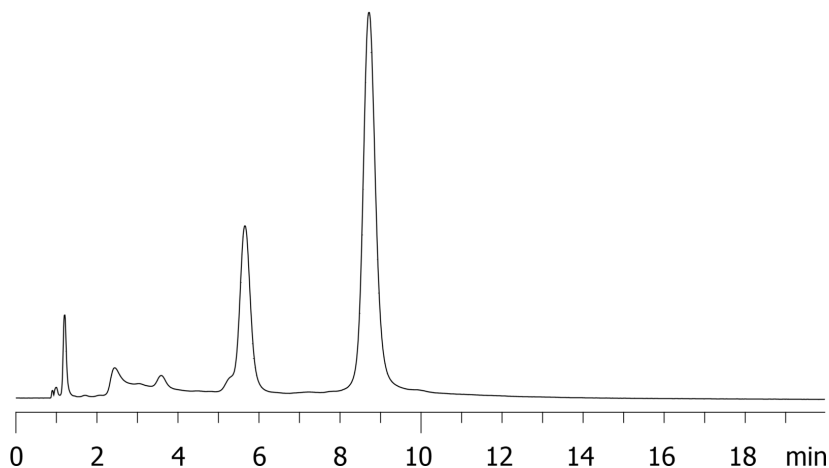
Gradient	Step No.	Time (min)	Pct A	Pct B	Pct C	Pct D
Profile:	1	0	70	10	20	0
	2	0	10	10	20	60

Flow Rate: 2.2 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Diode Array (PDA) @ (ambient)**Analyst Note:** SecurityGuardTM Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Clarity Oligo-WAX 4 x 3.0mm ID, 10/Pk Part No.: AJ0-8324

- Holder Part No.: KJ0-4282

16397

**ANALYTES:****1** AP site & DNA 20mer