

Oligo WAX - AP-A&G sites & DNA 20mer

Column: N/A
Dimensions: N/A
Order No: N/A
Elution Type: N/A
Eluent A: Water
Eluent B: Acetonitrile
Eluent C: 100mM Tris pH 8
Eluent D: 2M NaCl

Gradient	Step No.	Time (min)	Pct A	Pct B	Pct C	Pct D
Profile:	1	0	70	10	20	0
	2	20	10	10	20	60

Flow Rate: 2.2 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Diode Array (PDA) @ (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Clarity Oligo-WAX 4 x 3.0mm ID, 10/Pk Part No.: AJ0-8324

- Holder Part No.: KJ0-4282

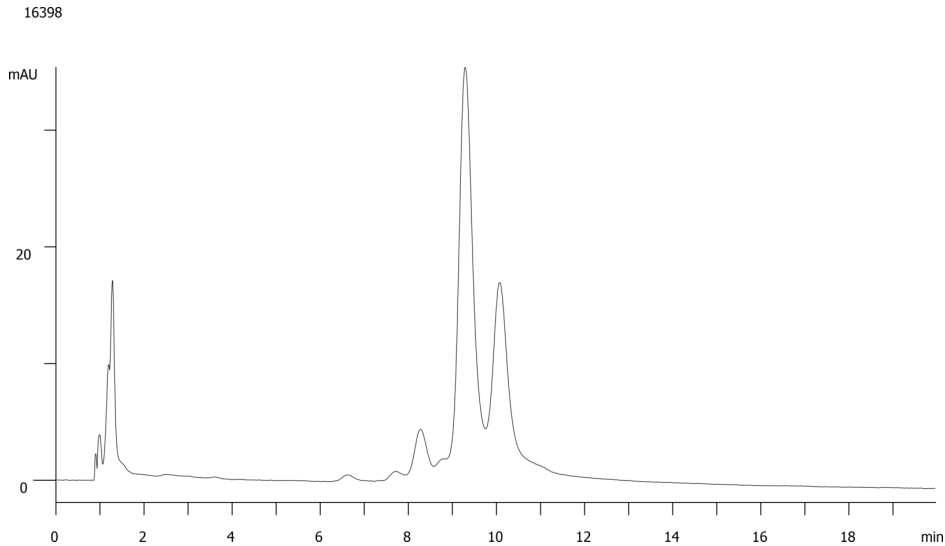

ANALYTES:

- 1 Depurinated A&G & 20mer DNA

