


EPA Method 552.2: Haloacetic Acids & Dalapon on ZB-XLB**Column:** Zebron™ ZB-XLB, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** Proprietary XLB Phase**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G019-11**Oven Profile:** 40 °C for 4min to 200 °C at 14 °C/min hold for 10min**Carrier Gas:** Constant Flow Helium, 1.3 mL/min**Injection:** Split 75:1 1 µL @ 200°C**Detection:** Electron Capture (ECD) (310°C)**Products used in this application:**

EPA Method 552.2: Haloacetic Acids & Dalapon on ZB-XLB

ANALYTES:

- 1 Chloroacetic acid
- 2 Bromoacetic acid
- 3 Dichloroacetic acid
- 4 Dalapon
- 5 Trichloroacetic acid
- 6 Bromochloroacetic acid
- 7 1,2,3-Trichloropropane
- 8 Bromodichloroacetic acid
- 9 Dibromoacetic acid
- 10 Chlorodibromoacetic acid
- 11 2,3-Dibromopropionic acid
- 12 Tribromoacetic acid

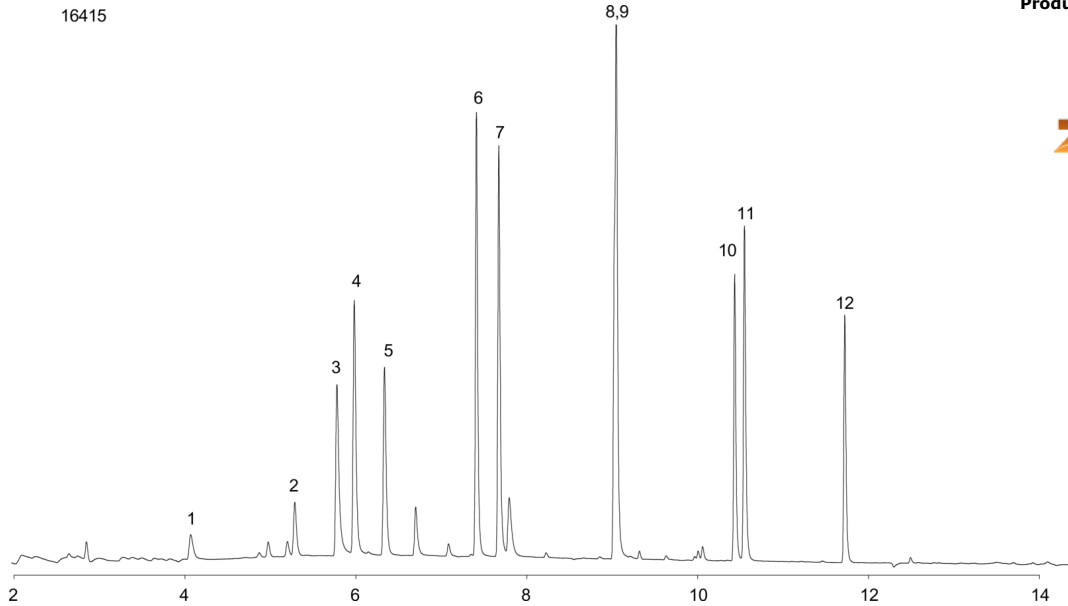