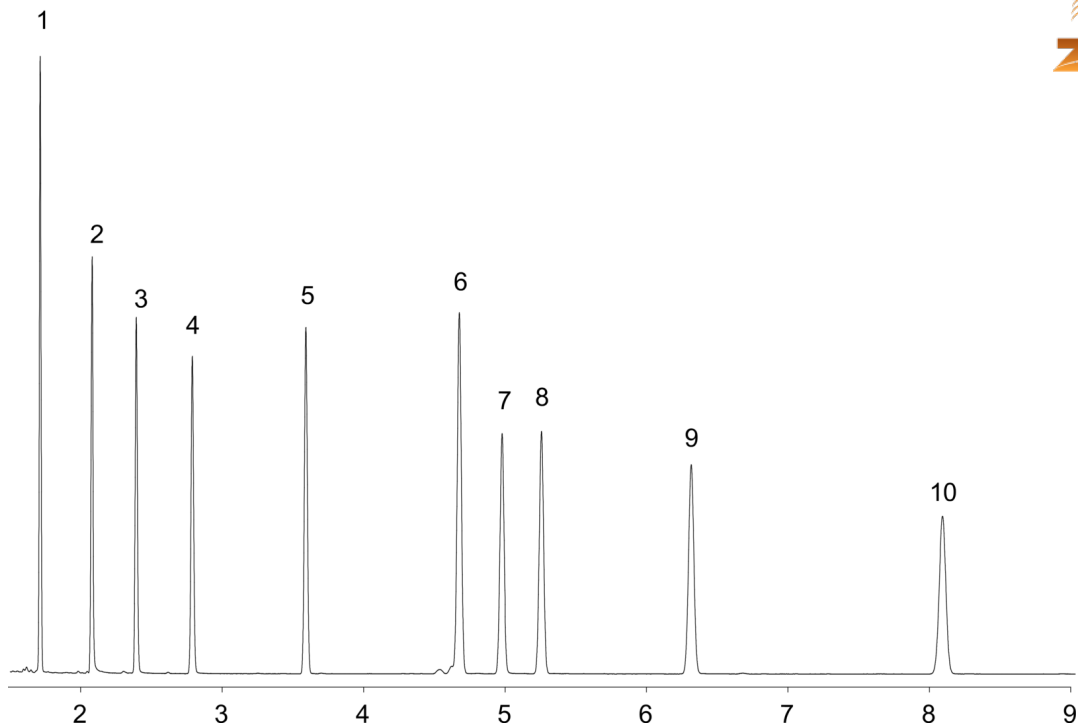


Industrial Chemicals Mix on ZB-XLB**Column:** Zebtron[™] ZB-XLB, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** Proprietary XLB Phase**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G019-11**Oven Profile:** 140 °C Isothermal**Carrier Gas:** Constant Flow Hydrogen, 1.2 mL/min**Injection:** Split 83:1 1 μ L @ 250°C**Detection:** Flame Ionization (FID) (325°C)

16417

**Products used in this application:**

Industrial Chemicals Mix on ZB-XLB

ANALYTES:

- 1 Decane
- 2 2-Ethylhexanoic acid
- 3 1,6-Hexanediol
- 4 4-Chlorophenol
- 5 Tridecane
- 6 1-Methylnaphthalene
- 7 1-Undecanol
- 8 Tetradecane
- 9 Dicyclohexylamine
- 10 Pentadecane

