


EPA 552.2 Haloacetic acids on ZB-XLB**Column:** Zebron[™] ZB-XLB, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** Proprietary XLB Phase**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G019-11**Oven Profile:** 48 °C for 15 min to 200 °C @ 10 °C/min.**Carrier Gas:** Constant Flow Helium, 1.5 mL/min**Injection:** Split 5:1 1 μ L @ 250°C**Detection:** Mass Selective (MSD) (180°C)**Analyst Note:** 50ppm methylated EPA 552.2 mix in methylene chloride.**Products used in this application:**

EPA 552.2 Haloacetic acids on ZB-XLB

ANALYTES:

- 1 Monochloroacetic acid
- 2 Monobromoacetic acid
- 3 Dichloroacetic acid
- 4 Dalapon
- 5 Trichloroacetic acid
- 6 Bromochloroacetic acid
- 7 1,2,3-Trichloropropane
- 8 Bromodichloroacetic acid
- 9 Dibromoacetic acid
- 10 Chlorodibromoacetic acid
- 11 2,3-Dibromopropionic acid
- 12 Tribromoacetic acid

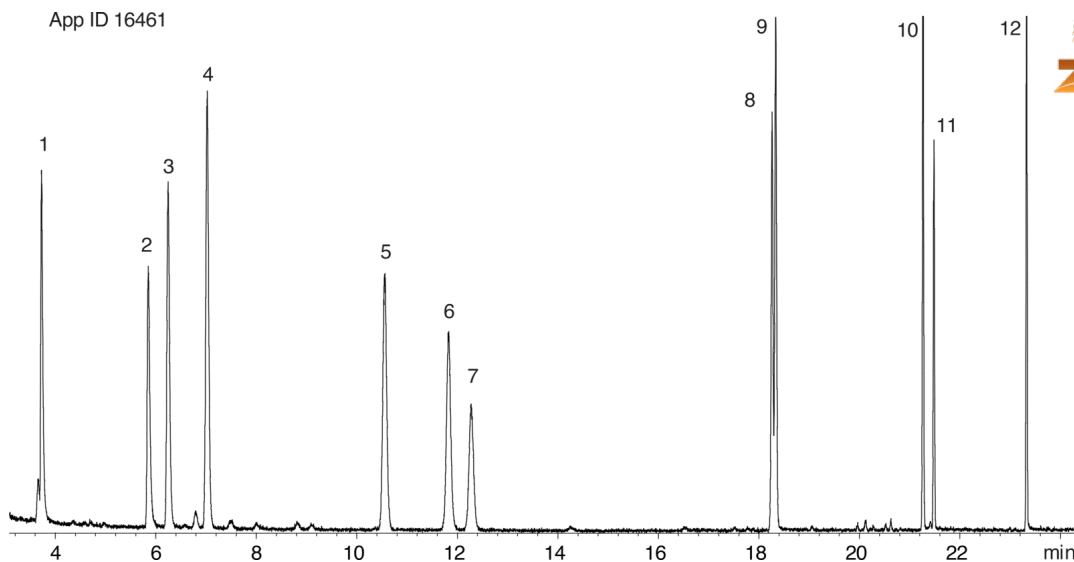