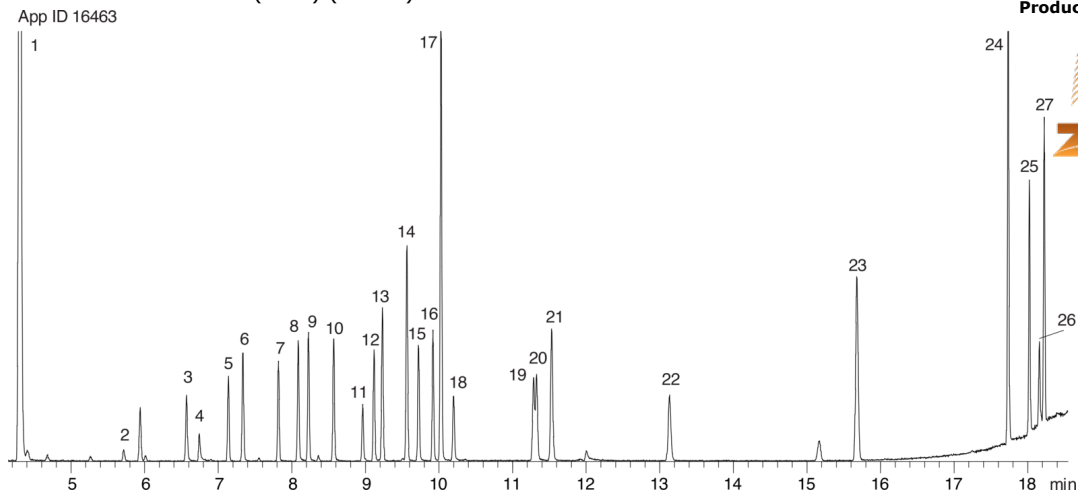


Drug Screen on ZB-XLB

Column: Zebron[™] ZB-XLB, GC Cap. Column 30 m x 0.25 mm x 0.50 µm, Ea
Phase: Proprietary XLB Phase
Dimensions: 30 meters x 0.25 mm x 0.5 µm
Order No: 7HG-G019-17
Oven Profile: 140 °C to 240°C @ 10°C/min for 5 min to 340 °C @ 25 °C/min for 5 min.

Carrier Gas: Constant Flow Helium, 1.5 mL/min
Injection: Split 5:1 1 µL @ 250°C
Detection: Mass Selective (MSD) (180°C)

**Products used in this application:**

Drug Screen on ZB-XLB

ANALYTES:

- 1 Nicotine
- 2 Methylecgonine
- 3 Acetaminophen
- 4 Ibuprofen
- 5 Allobarbitol
- 6 Aprobital
- 7 Butabital
- 8 Phenacetin
- 9 Amobarbitol
- 10 Pentobarbitol
- 11 Meprobamate
- 12 Secobarbitol
- 13 Methyl benzilate
- 14 Benzphetamine
- 15 Dimenhydrinate
- 16 Hexobarbitol
- 17 Caffeine
- 18 Doxylamine
- 19 Phenobarbitol
- 20 Methapyrilene
- 21 Chlorpheniramine
- 22 Brompheniramine
- 23 Cocaine
- 24 Chlorcyclizine
- 25 Codeine
- 26 Diazepam
- 27 Morphine
- 28 Hydrocodone

