


EPA 8151A Chlorinated Herbicides on ZB-XLB**Column:** Zebtron™ ZB-XLB, GC Cap. Column 30 m x 0.25 mm x 0.50 µm, Ea**Phase:** Proprietary XLB Phase**Dimensions:** 30 meters x 0.25 mm x 0.5 µm**Order No:** 7HG-G019-17**Oven Profile:** 60 °C to 300°C @ 10°C/min**Carrier Gas:** Constant Flow Helium, 1 mL/min**Injection:** Split 5:1 1 µL @ 240°C**Detection:** Mass Selective (MSD) (180°C)**Analyst Note:** 50ppm in methanol
IS/Surrogate approx. 10ppm**Products used in this application:**

EPA 8151A Chlorinated Herbicides on ZB-XLB

ANALYTES:

- 1 Dalapon
- 2 3,5-Dichlorobenzoic acid
- 3 4-Nitrophenol
- 4 DCAA
- 5 Dicamba
- 6 MCPP
- 7 MCPA
- 8 Dichlorprop
- 9 2,4-D
- 10 DBOB
- 11 Pentachlorophenol
- 12 Silvex
- 13 Chloramben
- 14 2,4,5-T
- 15 2,4-DB
- 16 Dinoseb
- 17 Bentazon
- 18 Picloram
- 19 Dacthal
- 20 Acifluorfen

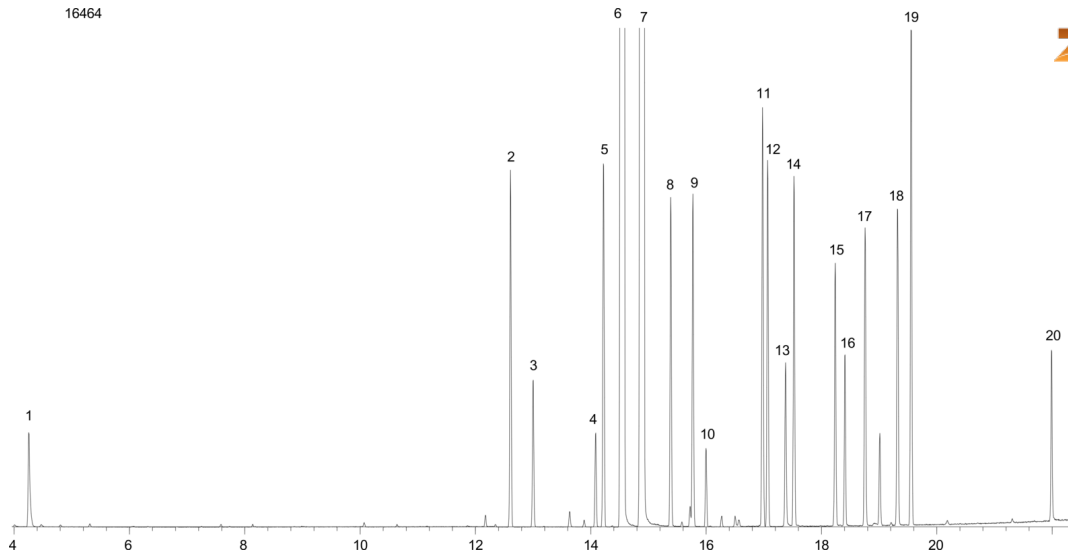