


Phenoxyacids Methyl Esters on ZB-MR1**Column:** Zebron[™] ZB-MultiResidue[™]-1, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G016-11**Oven Profile:** 60 °C for 0.5 min to 130 °C @ 10 °C/min to 210 °C @ 4 °C/min to 300 °C @ 10 °C/min**Carrier Gas:** Constant Flow Helium, 1.2 mL/min**Injection:** Splitless :1 1 μ L @ 260°C**Detection:** Mass Selective (MSD) (180°C)**Analyst Note:** Column with 5 meter guard. Analytes are 1.0 ppm.
Application courtesy of Annegret Tischendorf, Südsachsen Wasser GmbH**Products used in this application:**

Phenoxyacids Methyl Esters on ZB-MR1

ANALYTES:

- 1 Mecoprop
- 2 MCPA
- 3 Dichloroprop
- 4 2,4-D
- 5 Bromoxynil
- 6 Fenoprop
- 7 2,4,5-T
- 8 MCPB
- 9 2,4-DB
- 10 Bentazon

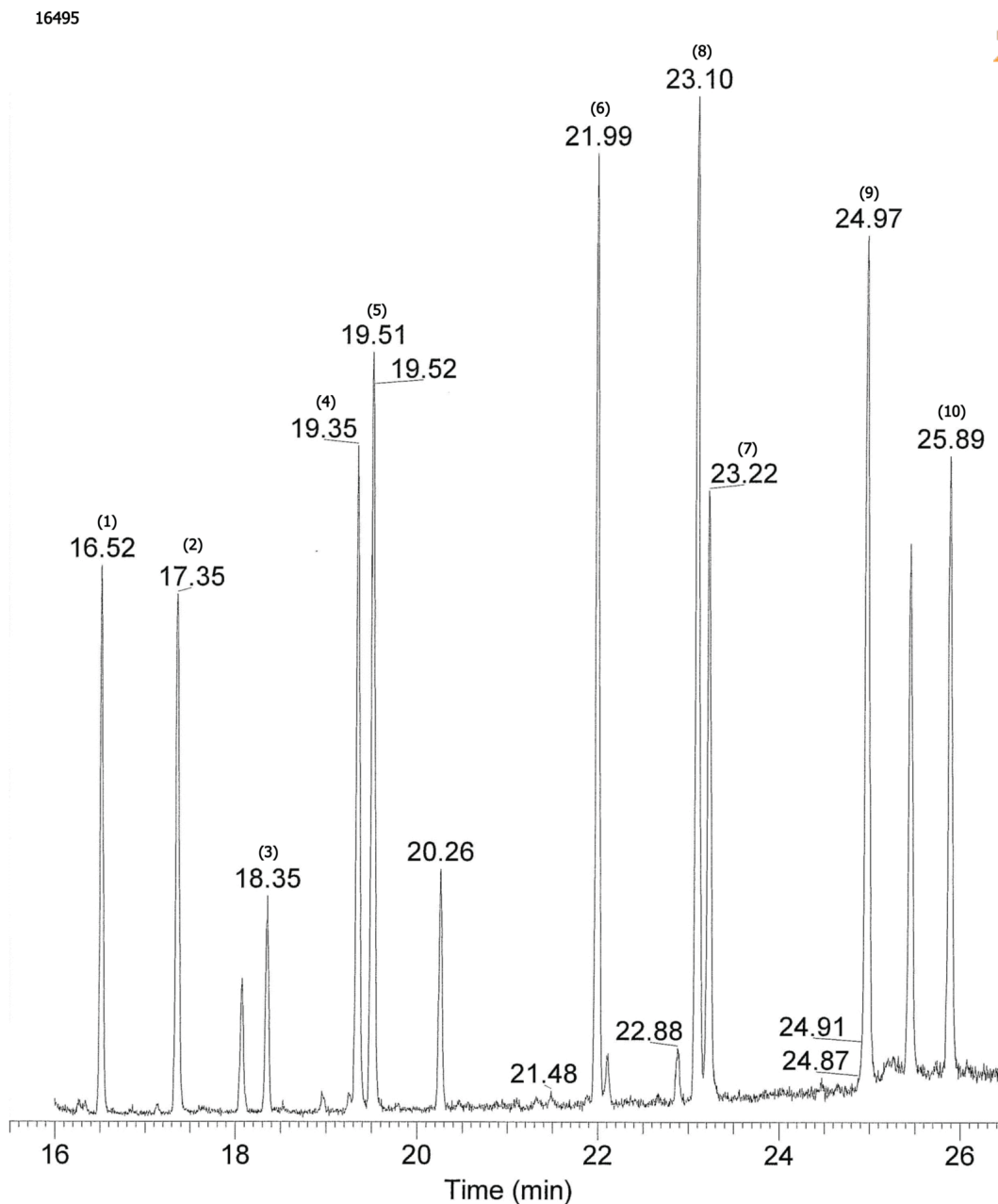