


TPH Diesel Analysis on ZB-1MS

Column: Zebron™ ZB-1MS, GC Cap. Column 20 m x 0.18 mm x 0.18 µm, Ea**Phase:** 100% Dimethylpolysiloxane**Dimensions:** 20 meters x 0.18 mm x 0.18 µm**Order No:** 7FD-G011-08**Oven Profile:** 75 °C for 1 min to 320 @ 30 °C/min for 3.83 min**Carrier Gas:** Ramp Pressure Helium, 2.5 mL/min**Injection:** Splitless :1 1 µL @ 280°C**Detection:** Flame Ionization (FID) (300°C)**Analyst Note:** Pressure ramp: 20. for 8 min to 3.0 mL/min @ 2.0 mL/min for 6 min.

Products used in this application:



16505

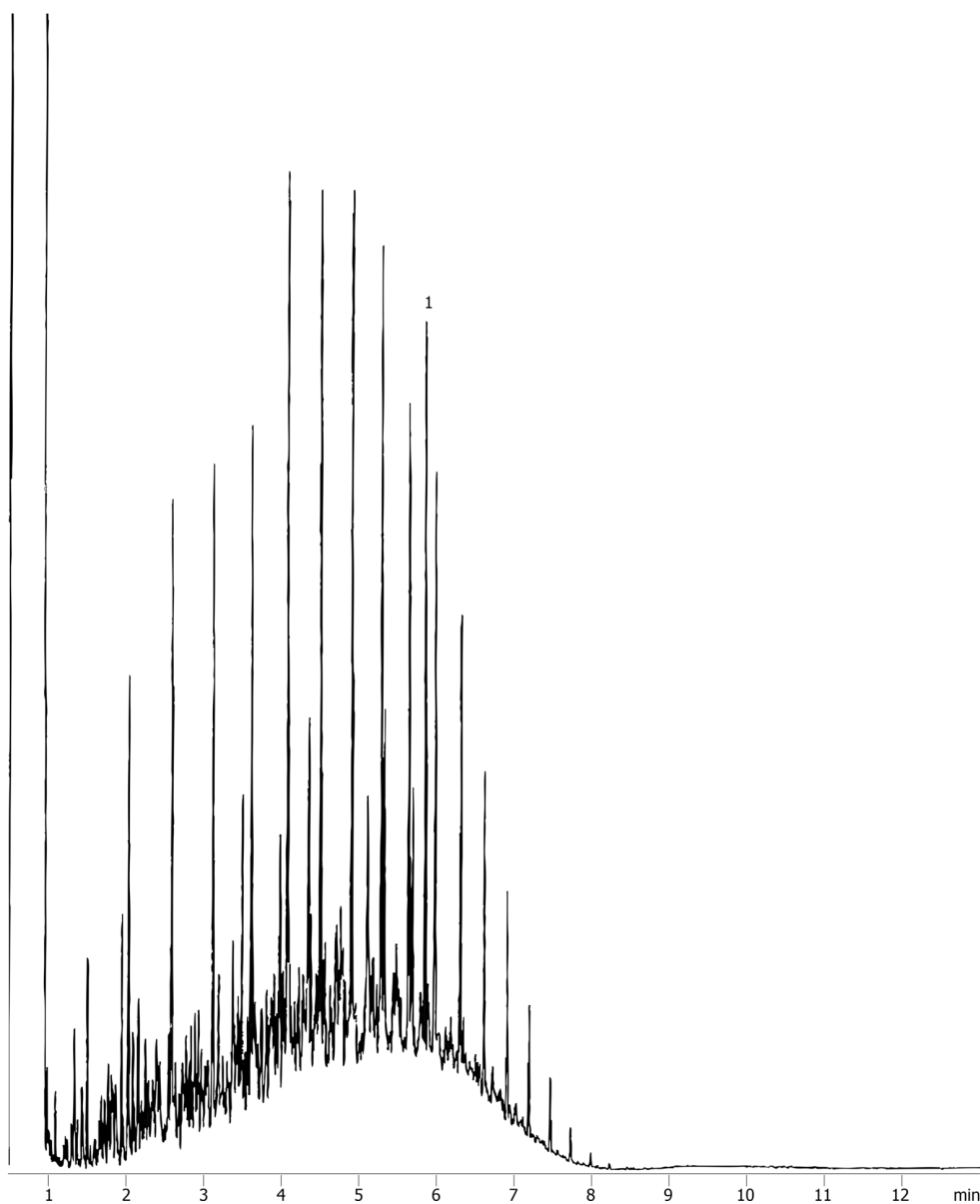

ANALYTES:

©2015 Phenomenex. All rights reserved.

For more information contact your Phenomenex Representative at info@phenomenex.com

Phenomenex products are available worldwide.

www.phenomenex.com.cninfo@phenomenex.com