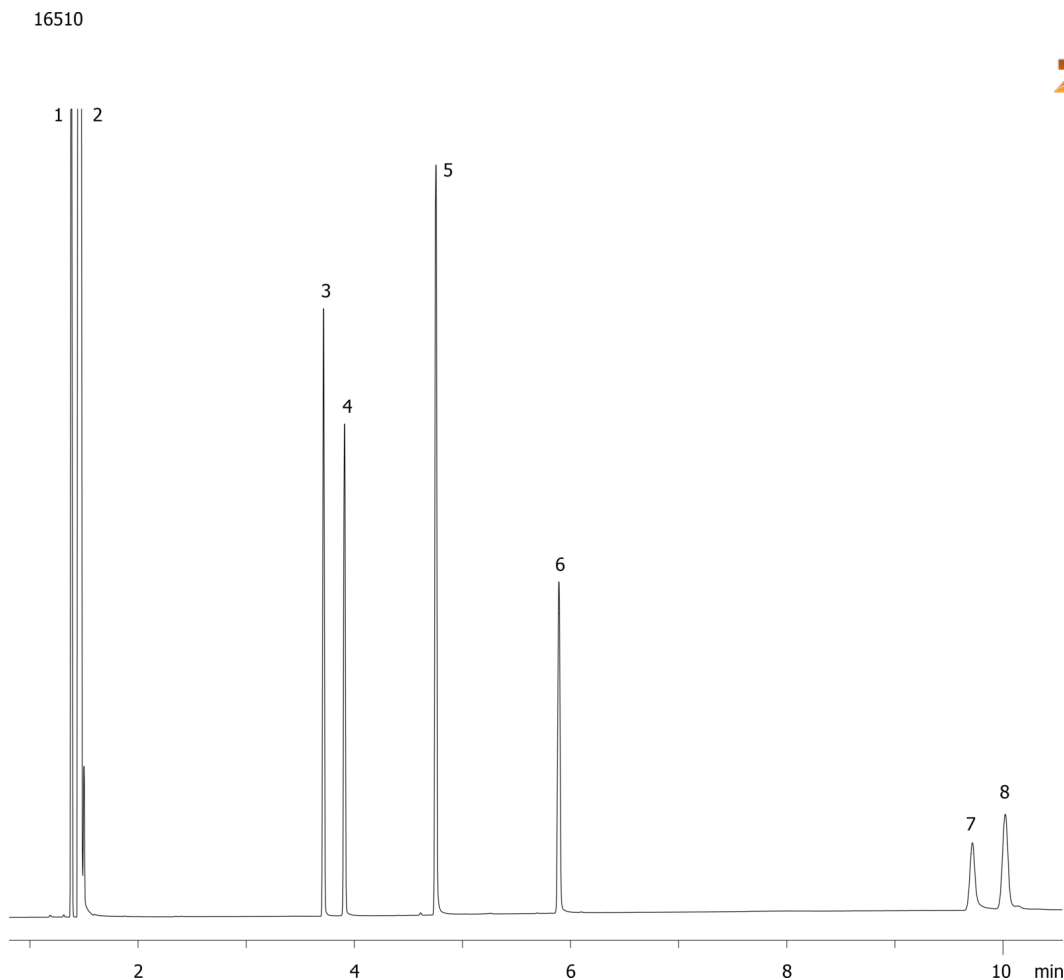


Glycols and Benzene in Toothpaste Standard on ZB-WAXplus**Column:** Zebron[™] ZB-WAXplus[™], GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** Polyethylene Glycol (PEG)**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G013-11**Oven Profile:** 80 °C to 180 °C @ 20 °C/min for 20 min.**Carrier Gas:** Constant Flow Helium, 3 mL/min**Injection:** Split 20:1 0.5 μ L @ 225°C**Detection:** Flame Ionization (FID) (250°C)**Analyst Note:** Analytes at 0.5 mg/ml**Products used in this application:****ANALYTES:**

- 1 Benzene
- 2 Acetonitrile
- 3 1,2-Propanediol
- 4 Ethylene glycol
- 5 1,3-Propanediol
- 6 Diethylene glycol
- 7 Glycerol
- 8 Triethylene glycol

©2005 Phenomenex, Inc. All rights reserved.

For more information contact your Phenomenex Representative at info@phenomenex.com

Phenomenex products are available worldwide.

www.phenomenex.com.cninfo@phenomenex.com