


Benzodiazepines by MSD on a ZB-XLB**Column:** Zebtron™ ZB-XLB, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** Proprietary XLB Phase**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G019-11**Oven Profile:** 180 °C for 0.5 min to 360 °C @ 20 °C/min.**Carrier Gas:** Constant Flow Helium, 1.5 mL/min**Injection:** Splitless :1 1 µL @ 250°C**Detection:** Mass Selective (MSD) (180°C)**Analyst Note:** 5ppm Benzodiazepines in methanol**Products used in this application:**

Benzodiazepines by MSD on a ZB-XLB

ANALYTES:

- 1 Oxazepam
- 2 Lorazepam
- 3 Diazepam
- 4 Nordazepam
- 5 Midazolam
- 6 Temazepam
- 7 Flunitrazepam
- 8 Prazepam
- 9 7-Aminoflunitrazepam
- 10 Bromazepam
- 11 Lormetazepam
- 12 Flurazepam
- 13 Estazolam
- 14 Alprazolam
- 15 Triazolam

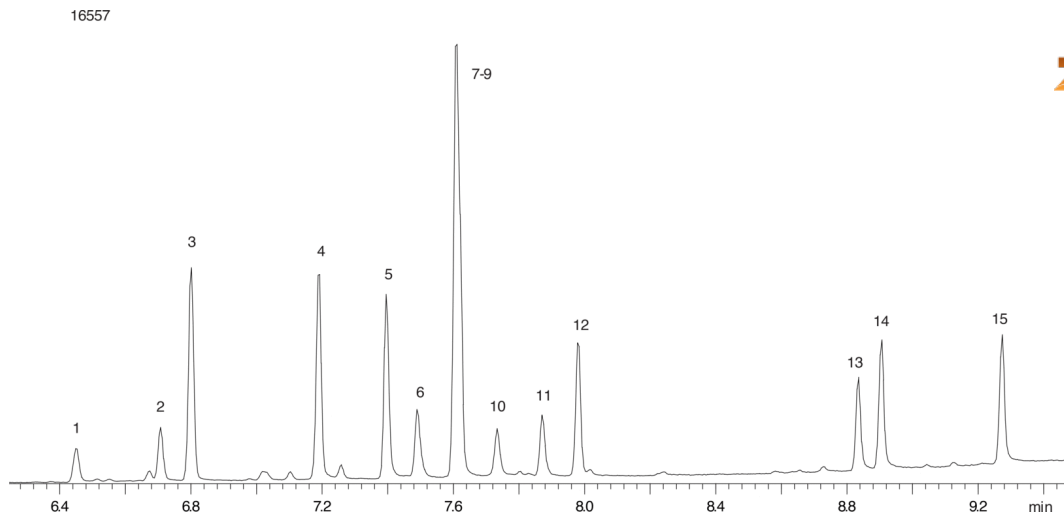