


EPA 8151A on MR2 (Carrier Gas: Hydrogen)**Column:** Zebtron[™] ZB-MultiResidue[™]-2, GC Cap. Column 30 m x 0.32 mm x 0.25 μ m, Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 30 meters x 0.32 mm x 0.25 μ m**Order No:** 7HM-G017-11**Oven Profile:** 50 °C for 0.5 min to 100 °C @ 25 °C/min to 205 °C at 22 °C/min to 320 °C @ 17 °C/min.**Carrier Gas:** Constant Flow Hydrogen, 2.4 mL/min**Injection:** Splitless :1 1 μ L @ 250°C**Detection:** Electron Capture (ECD) (335°C)**Analyst Note:** Analaytes are 0.5 ppm in Hexane**Products used in this application:**

EPA 8151A on MR2 (Carrier Gas: Hydrogen)

ANALYTES:

- 1 Dalapon
- 2 3,5-Dichlorobenzoic acid
- 3 4-Nitrophenol
- 4 DCAA
- 5 Dicamba
- 6 MCPP
- 7 MCPA
- 8 Dichloroprop
- 9 DBOB
- 10 2,4-D
- 11 contaminant
- 12 Pentachlorophenol
- 13 Silvex
- 14 Chloromben
- 15 2,4,5-T
- 16 Dinoseb
- 17 2,4-DB
- 18 Bentazon
- 19 Picloram
- 20 Dacthal
- 21 Acifluorfen

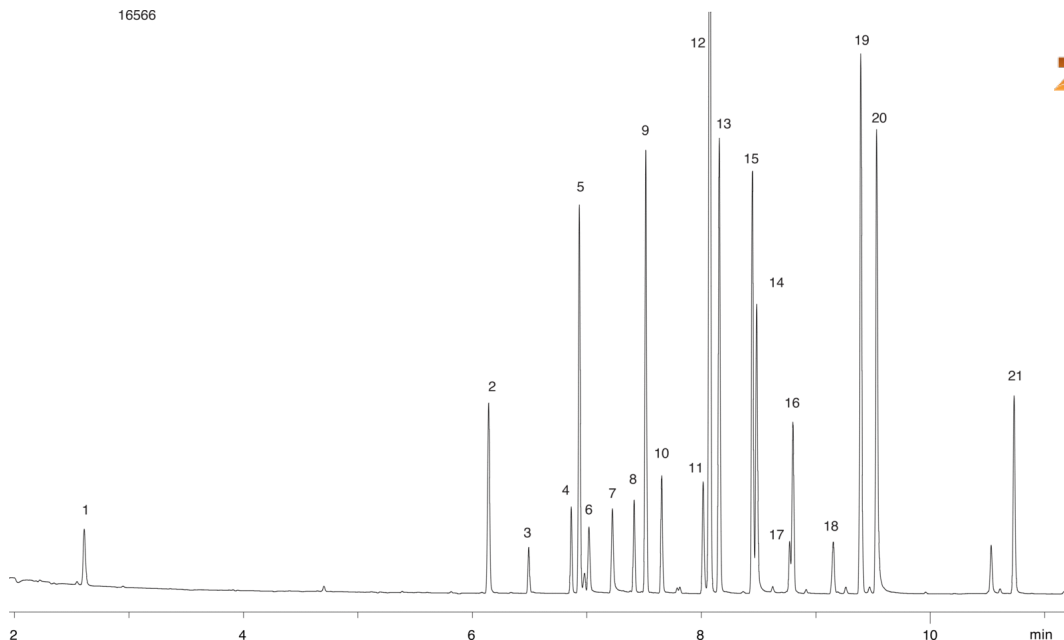