


EPA 552.2 on ZB-MR2 (Carrier Gas: Hydrogen)**Column:** Zebtron[™] ZB-MultiResidue[™]-2, GC Cap. Column 30 m x 0.32 mm x 0.25 μ m, Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 30 meters x 0.32 mm x 0.25 μ m**Order No:** 7HM-G017-11**Oven Profile:** 40 °C to 80 °C @ 12 °C/min for 4 min, to 225 °C @ 30 °C/min.**Carrier Gas:** Constant Flow Hydrogen, 2 mL/min**Injection:** Splitless :1 1 μ L @ 250°C**Detection:** Electron Capture (ECD) (250°C)**Analyst Note:** Standard concentration mix 0.5 ppm**Products used in this application:**

EPA 552.2 on ZB-MR2 (Carrier Gas: Hydrogen)

ANALYTES:

- 1 Chloroacetic acid
- 2 Bromoacetic acid
- 3 Dichloroacetic acid
- 4 Dalapon
- 5 Trichloroacetic acid
- 6 Bromochloroacetic acid
- 7 1,2,3-Trichloropropane
- 8 Dibromoacetic acid
- 9 Bromodichloroacetic acid
- 10 Chlorodibromoacetic acid
- 11 2,3-Dibromopropionic acid
- 12 Tribromoacetic acid

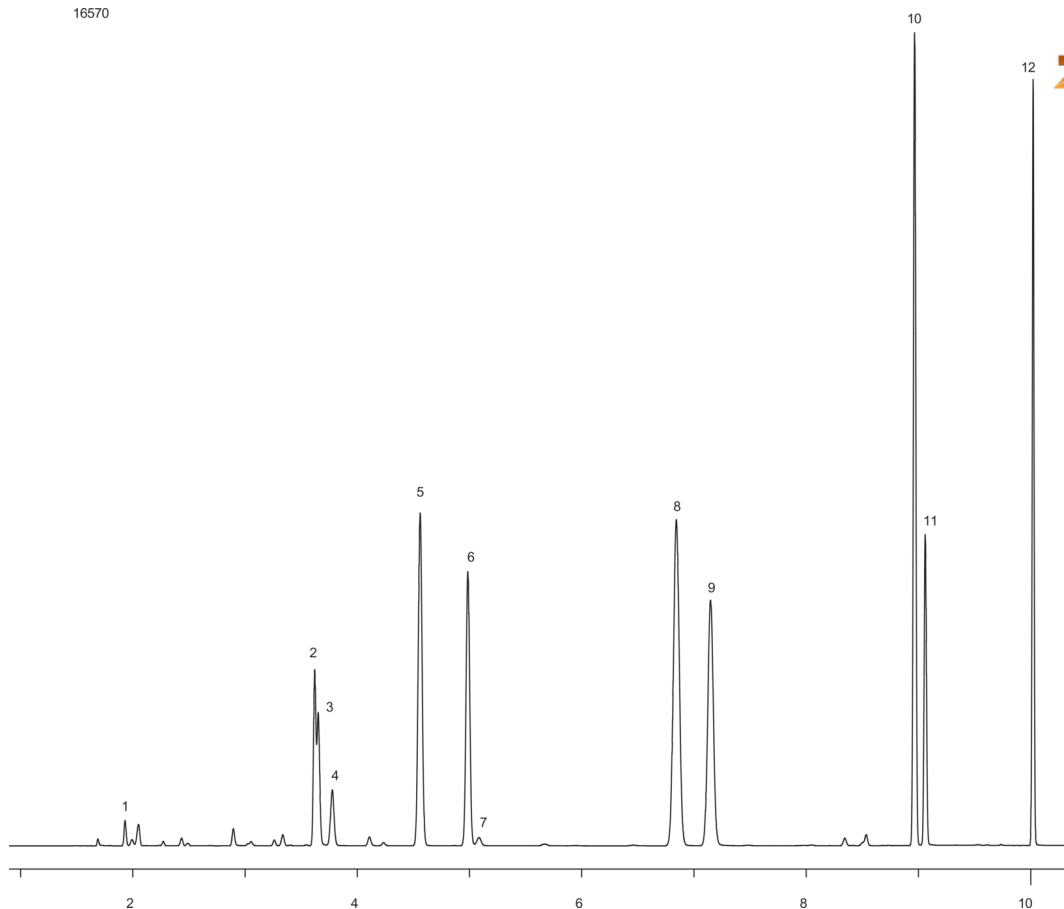