

Amitriptyline, Naproxen, Toluene in formic acid on Gemini-NX 5µm C18

Column: Gemini® 5 µm NX-C18 110 Å, LC Column 150 x 4.6 mm, Ea**Dimensions:** 150 x 4.6 mm ID**Order No:** 00F-4454-E0**Elution Type:** Gradient**Eluent A:** 0.1% Formic acid in water**Eluent B:** 0.1% Formic acid in ACN

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	95	5
	2	10	5	95
	3	12	5	95

Flow Rate: 1.5 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287


- Holder Part No.: KJ0-4282



Products used in this application:



16593

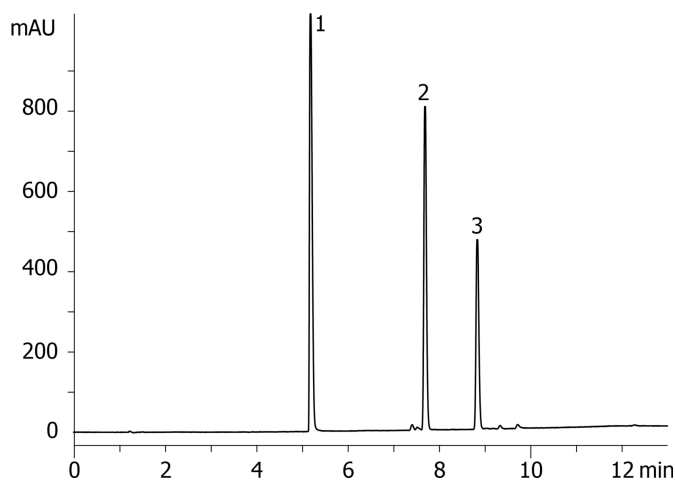

ANALYTES:

- 1 Amitriptyline
- 2 Naproxen
- 3 Toluene

