

Hydrophobic Bases (Diltiazem, Promethazine) in Formic Acid on Gemini-NX 5um C18

Column: Gemini® 5 µm NX-C18 110 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4454-E0

Elution Type: Gradient

Eluent A: 0.1% formic acid in water

Eluent B: 0.1% formic acid in acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	10	5	95

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)


Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282



Products used in this application:

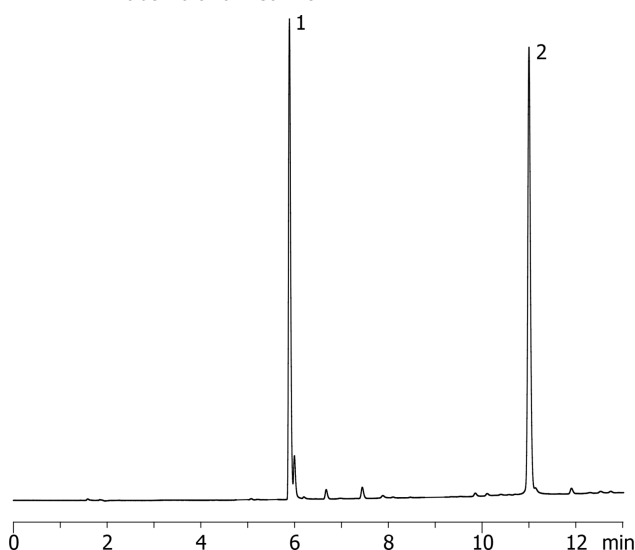

ANALYTES:

1 Diltiazem

Retention Time: 5.889 min

2 Promethazine

Retention Time: 11 min

