


ASTM D5501: Determination of Denatured Bioethanol on ZB-1**Column:** Zebron[™] ZB-1, GC Cap. Column 100 m x 0.25 mm x 0.50 μ m, Ea**Phase:** 100% Dimethylpolysiloxane**Dimensions:** 100 meters x 0.25 mm x 0.5 μ m**Order No:** 7MG-G001-17**Oven Profile:** 45 °C for 7 to 255 °C @ 30 °C/min for 6 min.**Carrier Gas:** Constant Flow Hydrogen, 35 cm/sec**Injection:** Split 50:1 1 μ L @ 300°C**Detection:** Flame Ionization (FID) (300°C)**Analyst Note:** Instrumentation: Shimadzu GC-2010

16664

**Products used in this application:****ANALYTES:**

- 1 Methanol
- 2 Ethanol

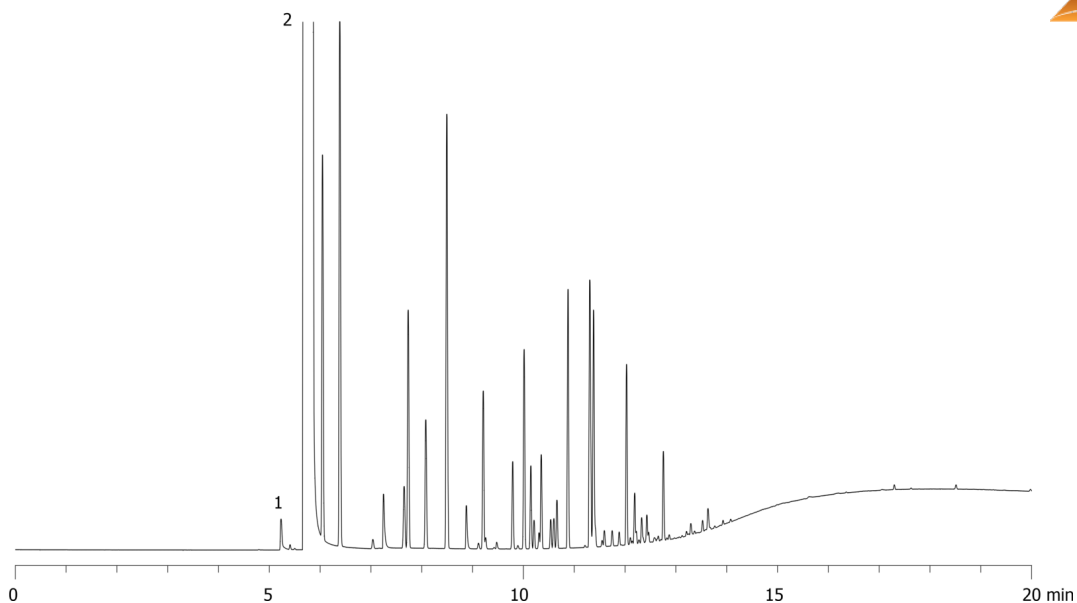