


ASTM D5501: Determination of Denatured Bioethanol on ZB-Bioethanol**Column:** Zebron[™] ZB-Bioethanol, GC Cap. Column 30 m x 0.25 mm x 1.00 μ m, Ea**Phase:** Unknown or Proprietary**Dimensions:** 30 meters x 0.25 mm x 1 μ m**Order No:** 7HG-G020-22**Oven Profile:** 45 °C for 2.5 min to 255 °C @ 30 °C/min for 4.5 min.**Carrier Gas:** Constant Flow Hydrogen, 35 cm/sec**Injection:** Split 50:1 1 μ L @ 300°C**Detection:** Flame Ionization (FID) (300°C)**Analyst Note:** Instrumentation: Shimadzu GC-2010**Products used in this application:**

16666

**ANALYTES:**

- 1 Methanol
- 2 Ethanol

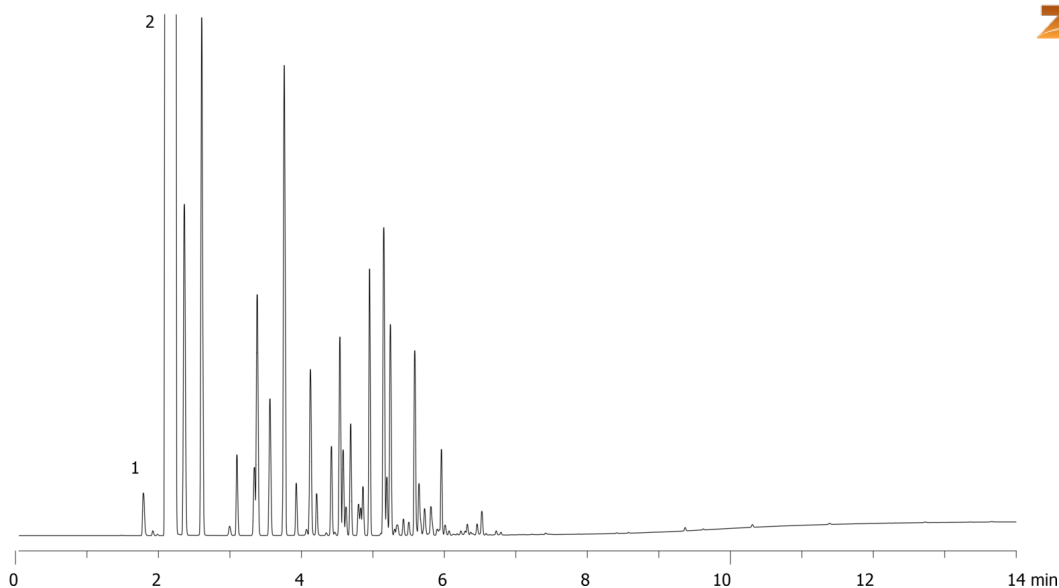