


Propoxyphene by MSD on a ZB-1MS**Column:** Zebtron[™] ZB-1MS, GC Cap. Column 10 m x 0.18 mm x 0.18 μ m, Ea**Phase:** 100% Dimethylpolysiloxane**Dimensions:** 10 meters x 0.18 mm x 0.18 μ m**Order No:** 7CD-G011-08**Oven Profile:** 150 °C to 320 °C @ 20 °C/min for 5 min**Carrier Gas:** Constant Flow Helium, 60 cm/sec**Injection:** Split 10:1 1 μ L @ 250°C**Detection:** Mass Selective (MSD) (320°C)**Analyst Note:** Analytes are 24ppm
Sim ions:
(Quant ion*) D11-Propoxyphene (*64,198,213): Propoxyphene (*58,193,208)**Products used in this application:****ANALYTES:**

- 1 Propoxyphene-d11
- 2 Propoxyphene

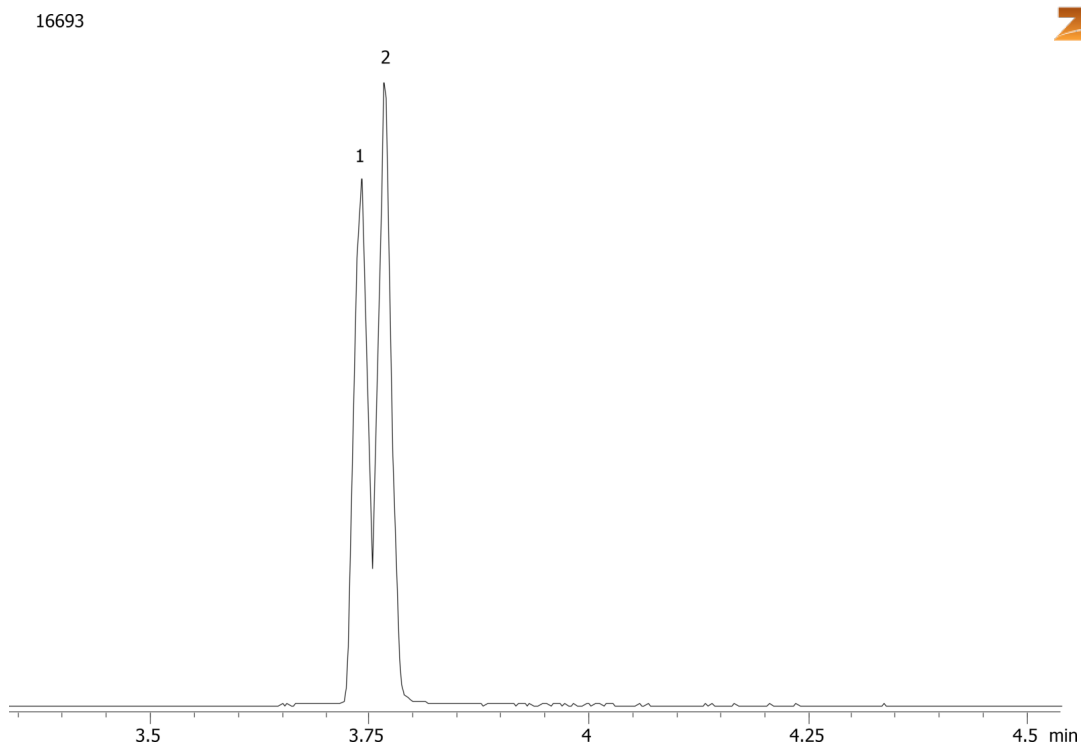