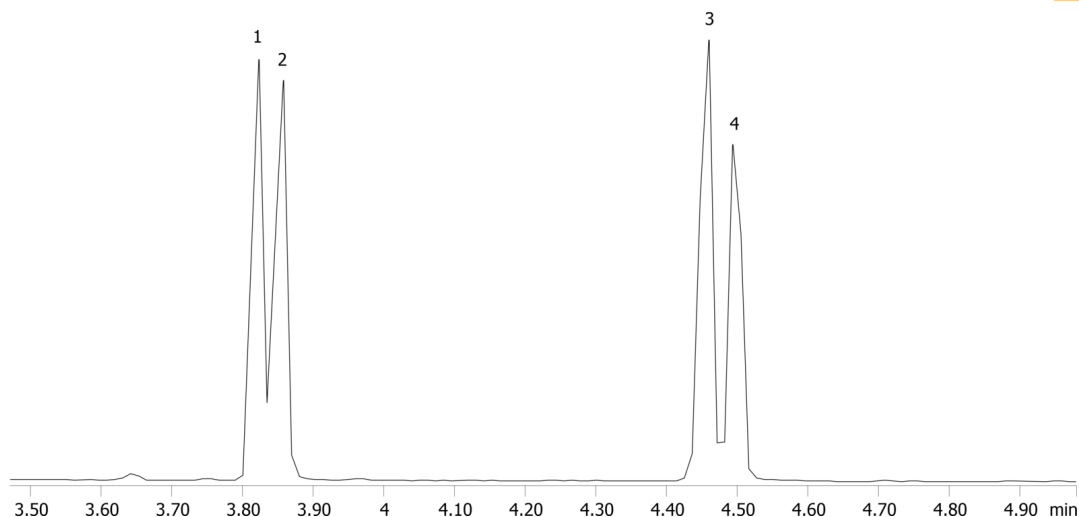


Amphetamines by MSD on a ZB-5MS**Column:** Zebron[™] ZB-5MS, GC Cap. Column 10 m x 0.18 mm x 0.18 µm, Ea**Phase:** 5% Phenyl-Arylene 95% Dimethylpolysiloxane**Dimensions:** 10 meters x 0.18 mm x 0.18 µm**Order No:** 7CD-G010-08**Oven Profile:** 50 °C to 320 °C @ 20 °C/ min**Carrier Gas:** Constant Flow He, 1.2 mL/min**Injection:** Split 10:1 1 µL @ 250°C**Detection:** Mass Selective (MSD) (320°C)**Analyst Note:** Analytes are 40ppm derivatized with HFBA

Sim ions:

(Quant ion*) D11-Amphetamine (*244,128): Amphetamine (*240,118,117): D14-Methamphetamine (*261,213):
Methamphetamine (*254,210,118)**Products used in this application:**

16737

**ANALYTES:**

- 1 Amphetamine-d11
- 2 Amphetamine
- 3 Methamphetamine-d14
- 4 Methamphetamine

