


Opiates by MSD on a ZB-5MS**Column:** Zebtron™ ZB-5MS, GC Cap. Column 10 m x 0.18 mm x 0.18 µm, Ea**Phase:** 5% Phenyl-Arylene 95% Dimethylpolysiloxane**Dimensions:** 10 meters x 0.18 mm x 0.18 µm**Order No:** 7CD-G010-08**Oven Profile:** 140 °C to 320 °C @ 20 °C/ min for 2 min.**Carrier Gas:** Constant Flow He, 1.2 mL/min**Injection:** Split 10:1 1 µL @ 250°C**Detection:** Mass Selective (MSD) (200°C)**Analyst Note:** Analytes are 10ppm derivatized with BSTFA
Sim ions:
(Quant ion*) D6-Codeine (*377,349): Codeine (*371,343,234): D3-Morphine (*432,417): Morphine (*429,414,401)**Products used in this application:**

16741

**ANALYTES:**

- 1 Codeine-d6
- 2 Codeine
- 3 Morphine-d3
- 4 Morphine

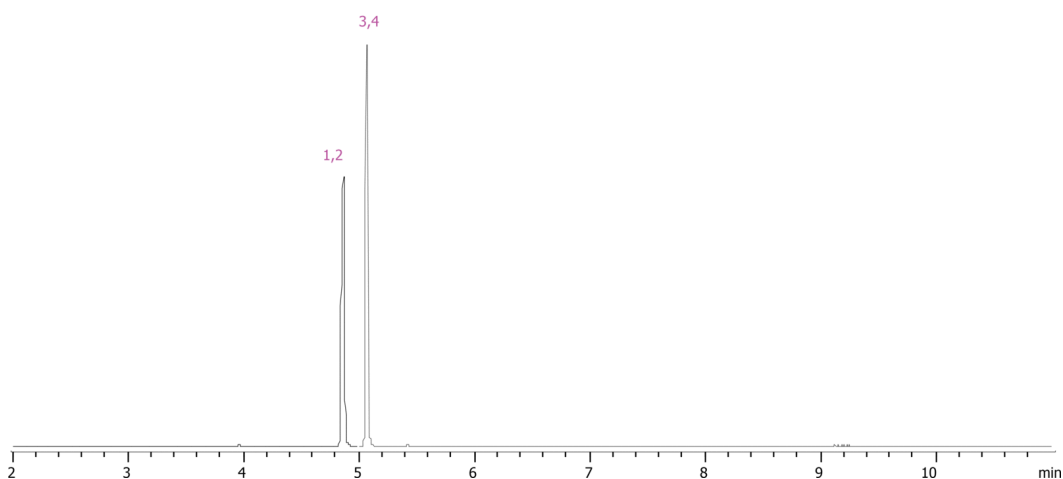