


**Sterols by FID on ZB-5HT****Column:** Zebtron™ ZB-5HT, GC Cap. Column 30 m x 0.25 mm x 0.10 µm, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.1 µm**Order No:** 7HG-G015-02**Oven Profile:** 220 °C to 350 °C @ 15 °C/min.**Carrier Gas:** Constant Flow Helium, 2 mL/min**Injection:** Splitless :1 0.5 µL @ 350°C**Detection:** Flame Ionization (FID) (350°C)**Analyst Note:** Sterols derivatized via BSTFA:TMCS; 99:1 in pyridine**Products used in this application:**

16749

**ANALYTES:**

- 1 Cholesterol
- 2 Campesterol
- 3 Stigmasterol
- 4 b-Sitosterol
- 5 Betulin

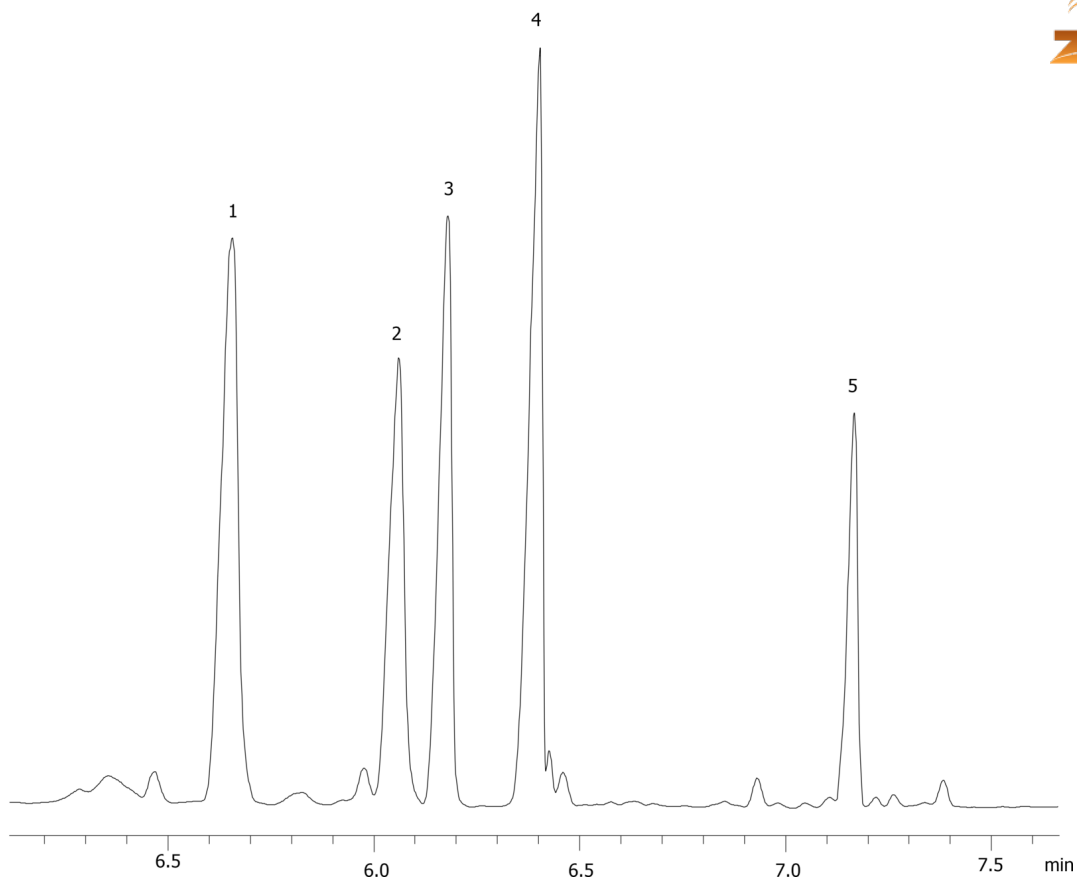