


Sterols from Peanut Butter by FID on ZB-MR2**Column:** Zebron™ ZB-MultiResidue™-2, GC Cap. Column 30 m x 0.25 mm x 0.20 µm, Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 30 meters x 0.25 mm x 0.2 µm**Order No:** 7HG-G017-10**Oven Profile:** 220 °C to 320 °C @ 13 °C/min.**Carrier Gas:** Constant Flow Helium, 2 mL/min**Injection:** Splitless :1 0.5 µL @ 300°C**Detection:** Flame Ionization (FID) (300°C)**Analyst Note:** Sample Prep: Saponification, Solid Phase Extraction (SPE), and Derivatization.**Products used in this application:****ANALYTES:**

- 1 Campesterol
- 2 Stigmasterol
- 3 b-Sitosterol
- 4 Betulin

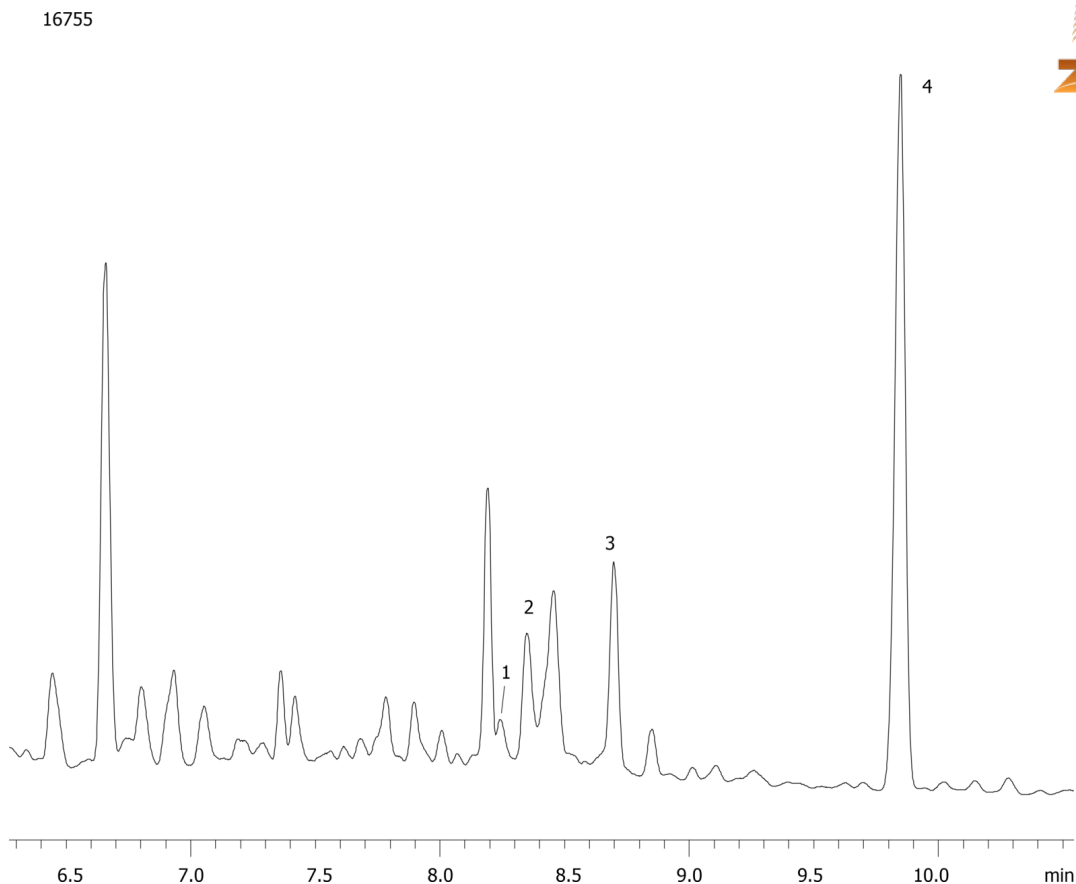