

Polar-RP - Analgesics

Column: Synergi™ 4 µm Polar-RP 80 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4336-E0

Elution Type: Isocratic

Eluent A: 20mM Ammonium formate pH 3.5

Eluent B: Acetonitrile


Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	50	50

Flow Rate: 2 mL/min

Col. Temp.: 25 °C

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 230 nm (ambient)

16899

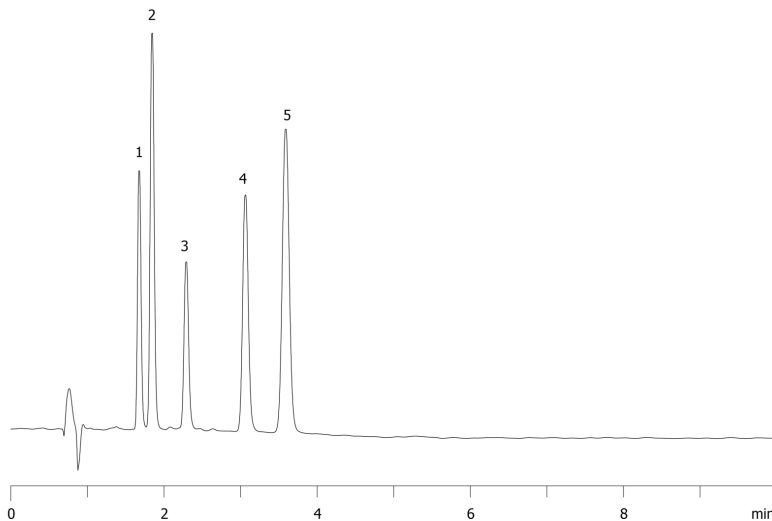

ANALYTES:

- 1 Ethylparaben
- 2 Indoprofen
- 3 Naproxen
- 4 Ibuprofen
- 5 Indomethacin

