

Warfarin on Lux Cellulose-2

Column: Lux® 5 µm Cellulose-2, LC Column 250 x 4.6 mm, Ea**Dimensions:** 250 x 4.6 mm ID**Order No:** 00G-4457-E0**Elution Type:** Isocratic**Eluent A:** 0.1% Formic acid in Hexane**Eluent B:** 0.1% Formic acid in Isopropanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	20	60	40

Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.


- SecurityGuard Cartridges, Lux Cellulose-2, 4 x 3.0mm, 10/Pk Part No.: AJ0-8366

- Holder Part No.: KJ0-4282

17312



Products used in this application:

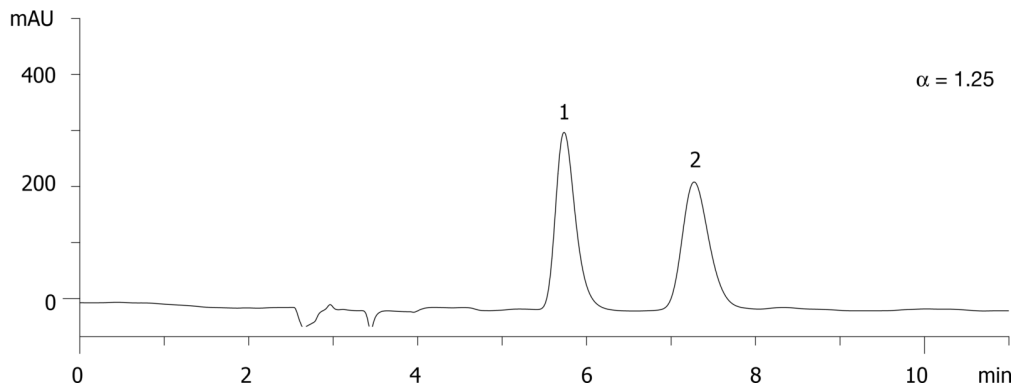

ANALYTES:

1 Warfarin

Retention Time: 5.7 min

2 Warfarin

Retention Time: 7.3 min

