

N-(3,5-Dinitrobenzoyl)-alpha-methylbenzylamine on Lux Cellulose-1

Column: Lux® 5 µm Cellulose-1, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4459-E0

Elution Type: Isocratic

Eluent A: 0.1% DEA in Hexane

Eluent B: 0.1% DEA in Isopropanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
1	1	20	60	40

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.


- SecurityGuard Cartridges, Lux Cellulose-1, 4 x 3.0mm, 10/Pk Part No.: AJ0-8403

- Holder Part No.: KJ0-4282

17386



Products used in this application:

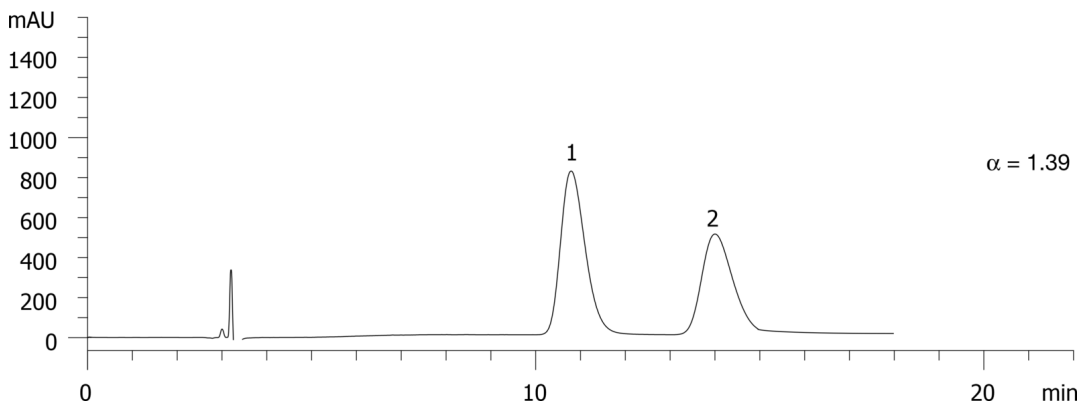

ANALYTES:

- 1 N-(3,5-Dinitrobenzoyl)-alpha-methylbenzylamine

Retention Time: 10.8 min

- 2 N-(3,5-Dinitrobenzoyl)-alpha-methylbenzylamine

Retention Time: 14 min

