

Lux Chiral Test Mix 5 (AL0-8412)

Column: Lux® 3 µm Cellulose-1, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4458-E0

Elution Type: Isocratic

Eluent A: Hexane

Eluent B: IPA

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
1	0	90	10	

Flow Rate: 0.5 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)

Analyst Note: This Check Standard is used to verify the performance of the following Lux chiral column phases: Cellulose 1 and 2. Order part number AL0-8412



Products used in this application:

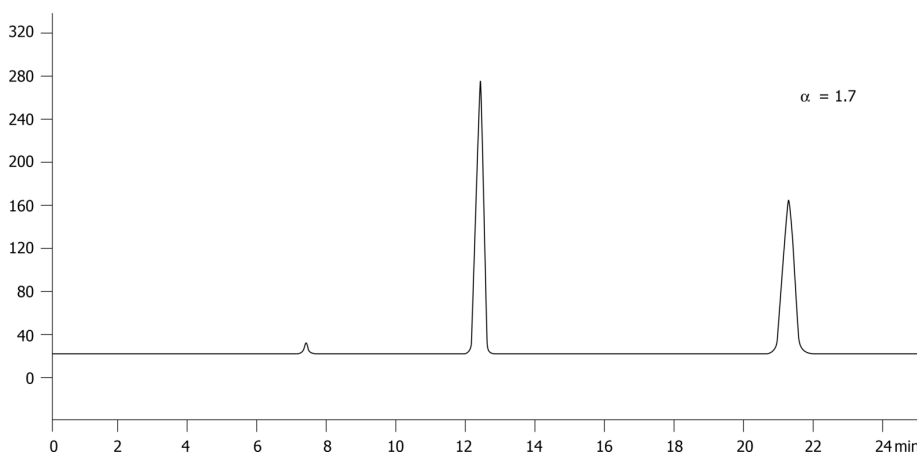

ANALYTES:

1 trans-Stilbene oxide

Retention Time: 12.35 min

2 trans-Stilbene oxide

Retention Time: 21.25 min

