

Milnacipran on Lux Cellulose-1**Column:** Lux® 5 µm Cellulose-1, LC Column 250 x 4.6 mm, Ea**Dimensions:** 250 x 4.6 mm ID**Order No:** 00G-4459-E0**Elution Type:** Isocratic**Eluent A:** 0.1 % DEA in Hexane**Eluent B:** 0.1% DEA in IPA

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	20	80	20

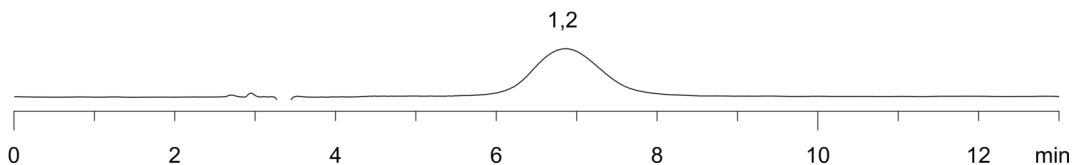
Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Cellulose-1, 4 x 3.0mm, 10/Pk Part No.: AJ0-8403

- Holder Part No.: KJ0-4282

**Products used in this application:**

17499

**ANALYTES:****1** Milnacipran

Retention Time: 6.8 min

