

Clenbuterol on Lux Cellulose-2**Column:** Lux® 5 µm Cellulose-2, LC Column 250 x 4.6 mm, Ea**Dimensions:** 250 x 4.6 mm ID**Order No:** 00G-4457-E0**Elution Type:** Isocratic**Eluent A:** 0.1 % DEA in CH₃CN**Eluent B:** 0.1% DEA in Water

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	20	60	40

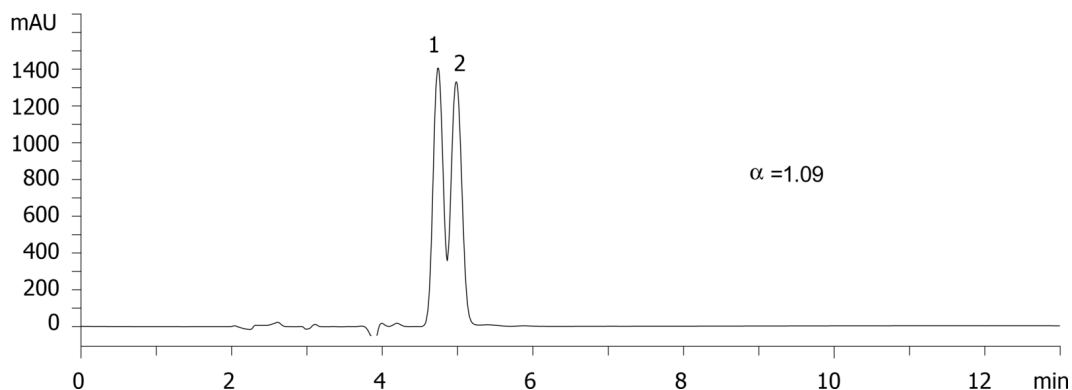
Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Cellulose-2, 4 x 3.0mm, 10/Pk Part No.: AJ0-8366

- Holder Part No.: KJ0-4282

**Products used in this application:**

17505

**ANALYTES:****1** Clenbuterol

Retention Time: 4.7 min

2 Clenbuterol

Retention Time: 4.9 min

