

## PEGylated proteins on Jupiter® 300 C18 3µm and 5µm, and C4 5µm columns

**Column:** Jupiter® 3 µm C18 300 Å, LC Column 150 x 4.6 mm, Ea

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-4263-E0

**Elution Type:** Gradient

**Eluent A:** 0.1% TFA, 2% ACN in Water

**Eluent B:** 90% ACN, 0.085% TFA in Water





| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B |
|-------------------|----------|------------|-------|-------|
|                   | 1        | 0          | 80    | 20    |
|                   | 2        | 25         | 35    | 65    |

**Flow Rate:** 1 mL/min

**Col. Temp.:** 45 °C

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 214 nm (25 °C)

**Analyst Note:** Compare reverse phase separations of different PEGylated and intact proteins selected at different reaction time-points on three different columns.

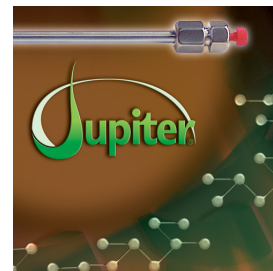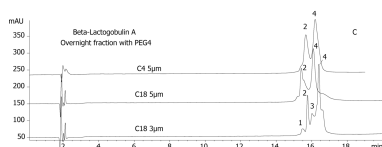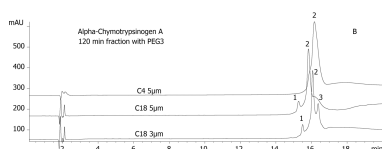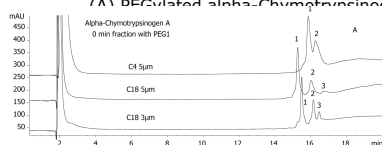

Products used in this application:



### ANALYTES:

- 1 PEGylated Proteins

