


**Organophosphorous Pesticides Mix 2 on ZB-MR1****Column:** Zebtron™ ZB-MultiResidue™-1, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G016-11**Oven Profile:** 40 °C for 1 min to 100 °C @ 15 °C/min to 225 °C @ 5 °C/min to 320 °C @ 15 °C/min for 5 min**Carrier Gas:** Constant Flow Helium, 0.7 mL/min**Injection:** On-Column :1 0.5 µL @ 43°C**Detection:** Mass Selective (MSD) (335°C)**Analyst Note:** Analytes 10 ppm in Acetone**Products used in this application:**

## Organophosphorous Pesticides Mix 2 on ZB-MR1

### ANALYTES:

- 1 Acephate
- 2 Omethoate
- 3 Cadusafos
- 4 Monocrotophos
- 5 Dimethoate
- 6 Tris(chloropropyl) phosphate
- 7 Tris(chloropropyl) phosphate
- 8 Tolclofos-methyl
- 9 Methyl parathion
- 10 Pirimiphos-methyl
- 11 Fenthion
- 12 Parathion
- 13 Quinalphos
- 14 Methidathion
- 15 Prothiophos
- 16 Tribufos
- 17 Triazophos
- 18 EPN
- 19 Phosalone

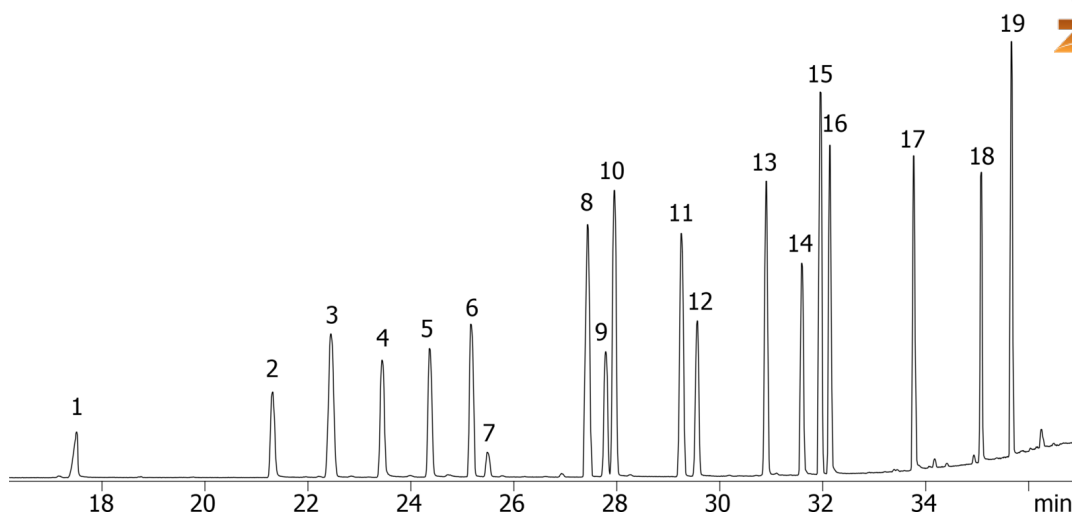