


Organophosphorous Pesticides Mix 3 on ZB-MR1**Column:** Zebron[™] ZB-MultiResidue[™]-1, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G016-11**Oven Profile:** 40 °C for 1 min to 100 °C @ 15 °C/min to 225 °C @ 5 °C/min to 320 °C @ 15 °C/min for 5 min**Carrier Gas:** Constant Flow Helium, 0.7 mL/min**Injection:** On-Column :1 0.5 μ L @ 43°C**Detection:** Mass Selective (MSD) (335°C)**Analyst Note:** Analytes 10 ppm in Acetone**Products used in this application:**

Organophosphorous Pesticides Mix 3 on ZB-MR1

ANALYTES:

- 1 Heptenophos
- 2 Ethoprop
- 3 Dyfonate
- 4 Cyanophos
- 5 Methyl paraoxon
- 6 Malaoxon
- 7 Oxydemeton methyl sulfone
- 8 Phorate sulfoxide
- 9 Phorate Sulfone
- 10 Isofenphos
- 11 Tetrachlorvinphos
- 12 Edifenphos
- 13 Fenamiphos
- 14 Fenamiphos sulfone
- 15 Bolstar

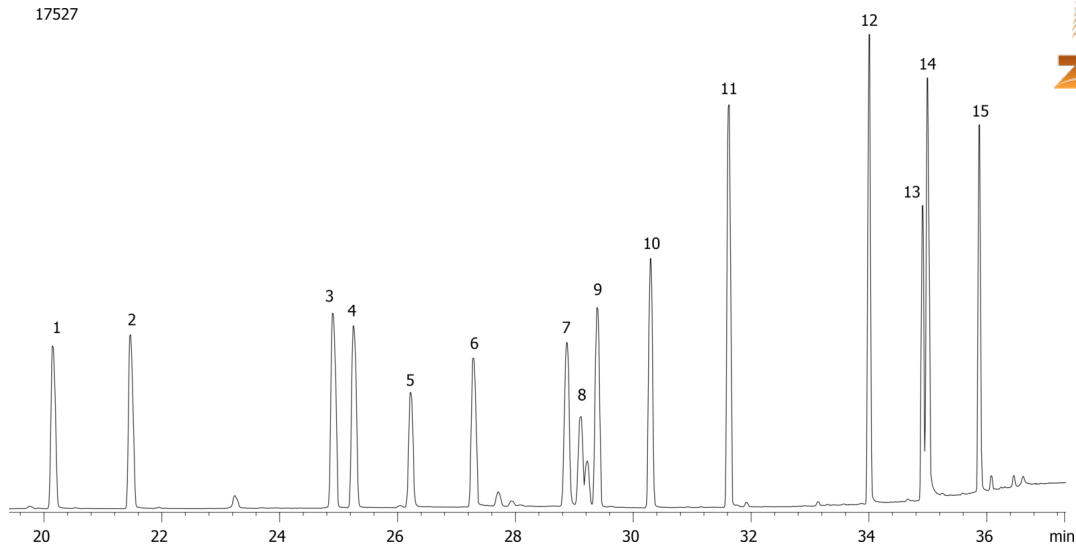