


NSO (Nitrogen/Sulfur/Oxygen) Pesticides Mix # 5 on ZB-MR1**Column:** Zebron[™] ZB-MultiResidue[™]-1, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G016-11**Oven Profile:** 40 °C for 1 min to 100 °C @ 15 °C/min to 225 °C @ 5 °C/min to 320 °C @ 15 °C/min for 5 min**Carrier Gas:** Constant Flow Helium, 0.7 mL/min**Injection:** On-Column :1 0.5 μ L @ 43°C**Detection:** Mass Selective (MSD) (335°C)**Analyst Note:** Analytes 10 ppm in Acetone**Products used in this application:**

17530

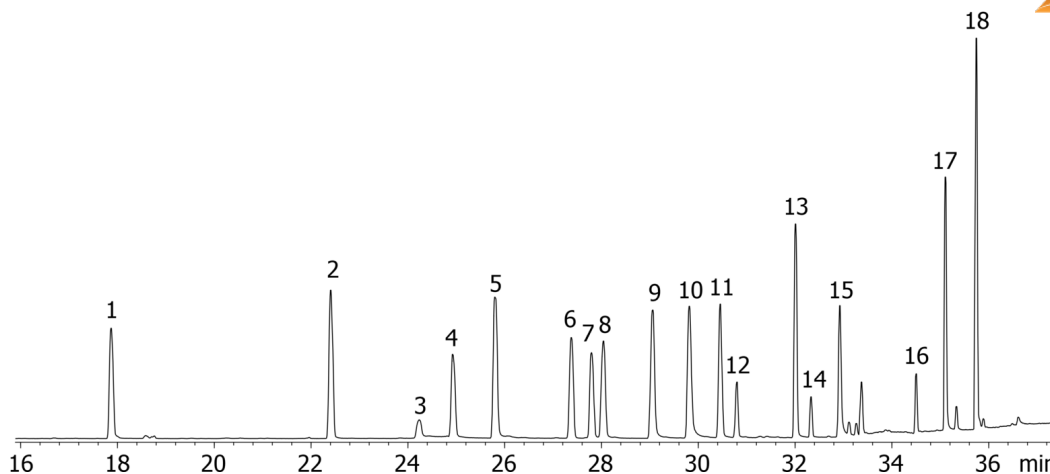

NSO (Nitrogen/Sulfur/Oxygen) Pesticides Mix # 5 on ZB-MR1

ANALYTES:

- 1 Metolcarb
- 2 2,3,5-Trimethacarb
- 3 Desmedipham
- 4 Cycluron
- 5 Secbumetone
- 6 Dioxacarb
- 7 Prometryne
- 8 3-Hydroxycarbofuran
- 9 Naphaleneacetamide
- 10 Tetraconazole
- 11 Pyracarbolid
- 12 Dinobuton
- 13 Napropamide
- 14 Pyrethrin 1
- 15 Tricyclazole
- 16 Nitralin
- 17 Fenoxycarb
- 18 Amitraz

