


Halogenated Pesticide Mix # 1 on ZB-MR1**Column:** Zebron[™] ZB-MultiResidue[™]-1, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G016-11**Oven Profile:** 40 °C for 1 min to 100 °C @ 15 °C/min to 225 °C @ 5 °C/min to 320 °C @ 15 °C/min for 5 min**Carrier Gas:** Constant Flow Helium, 0.7 mL/min**Injection:** On-Column :1 0.5 μ L @ 43°C**Detection:** Mass Selective (MSD) (335°C)**Analyst Note:** Analytes 10 ppm in Acetone**Products used in this application:**

Halogenated Pesticide Mix # 1 on ZB-MR1

ANALYTES:

- 1 Dichlobenil
- 2 BHC-gamma
- 3 Chlorothalonil
- 4 p,p'-Dicofol IP frag
- 5 Triadimenol
- 6 Captan
- 7 Folpet
- 8 Imazalil
- 9 p,p'-DDE
- 10 o,p'-Methoxychlor
- 11 Endosulfan sulfate
- 12 Methoxychlor
- 13 Fenhexamid
- 14 Dicofol
- 15 Acetamiprid
- 16 L-Cyhalothrin
- 17 Pyridaben
- 18 Cypermethrin
- 19 Cypermethrin
- 20 Cypermethrin
- 21 Pyraclostrobin

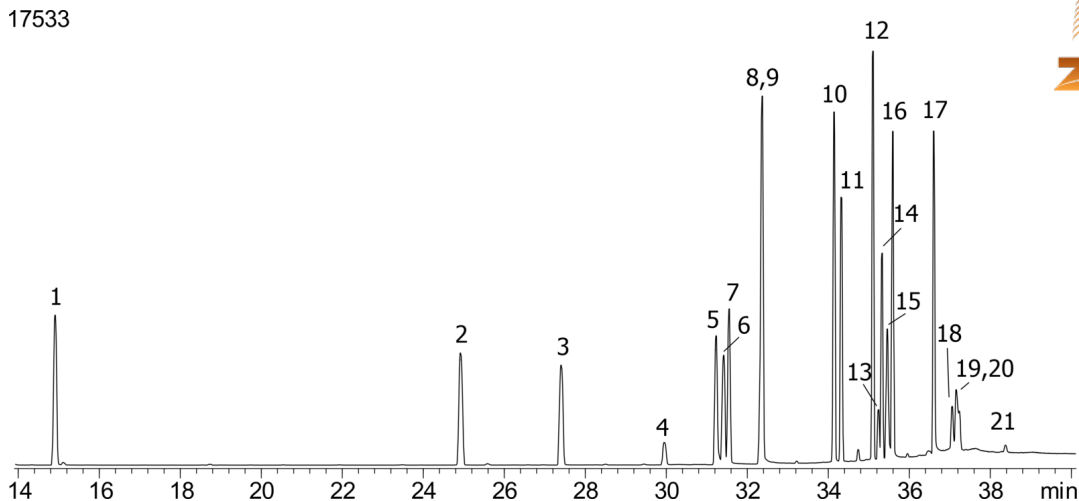